

Rice Lake TE Series

Tuning Fork Enhanced Balance

Operation Manual



© Rice Lake Weighing Systems. All rights reserved.

Rice Lake Weighing Systems® is a registered trademark of
Rice Lake Weighing Systems.

All other brand or product names within this publication are trademarks or
registered trademarks of their respective companies.

All information contained within this publication is, to the best of our knowledge, complete and
accurate at the time of publication. Rice Lake Weighing Systems reserves the right to make
changes to the technology, features, specifications and design of the equipment without notice.

The most current version of this publication, software, firmware and all other product
updates can be found on our website:

www.ricelake.com

Revision History

This section tracks and describes manual revisions for awareness of major updates.

Revision	Date	Description
D	August 22, 2025	Established revision history, general operation edits

Table i. Revision Letter History



Technical training seminars are available through Rice Lake Weighing Systems. Course descriptions and dates can be viewed at www.ricelake.com/training or obtained by calling 715-234-9171 and asking for the training department.

Contents

1.0	Introduction	7
1.1	Safety	8
1.2	FCC Compliance	8
1.3	Ensuring Accurate Measurements	9
1.3.1	Measuring Environment Precautions	9
1.3.2	Specimen Related Precautions	10
1.3.3	Balance Main Unit Precautions	10
2.0	Installation	11
2.1	Unpack Balance	11
2.1.1	Connections and Part Locations	12
2.2	Assemble Balance	13
2.2.1	Round Pan Base Assembly	13
2.2.2	Square Pan Base Assembly	13
2.3	Battery Installation	14
2.4	Leveling the Balance	15
2.5	Draft Shield Assembly	16
2.5.1	Draft Shield Assembly Parts	20
3.0	Operation	21
3.1	Front Panel	21
3.1.1	NTEP Only Models	21
3.1.2	NTEP / Measurement Canada Models	21
3.1.3	Front Panel Item Descriptions	22
3.2	Power On/Off	22
3.2.1	Standby	23
3.2.2	Balance Operation Check	23
3.3	Display Navigation	23
3.4	Zero-Point Adjustment	23
3.4.1	Zero-Point Adjustment Limits	23
3.5	Tare Value	24
3.5.1	Preset Tare	24
3.5.2	Check Tare Weight Using an F-Key	24
3.5.3	Add to Product	24
3.6	Weigh Mode	25
3.7	Counting Mode	25
3.7.1	Actual Value Setting	26
3.7.2	Numeric Value Setting Method	26
3.8	Percentage Mode	27
3.9	Coefficient Mode	28
3.10	Specific Gravity Mode	29
3.10.1	Menus Available in Operation Mode	29
3.10.2	Materials Required in Specific Gravity Mode	29
3.10.3	Prepare Equipment to Measure Specific Gravity	29
3.10.4	Measure the Specific Gravity Mode	29
3.11	Statistics Mode	30
3.11.1	Menus Available in Operation Mode	30



Rice Lake continually offers web-based video training on a growing selection of product-related topics at no cost. Visit www.ricelake.com/webinars

3.11.2	Selecting Statistic Mode	31
3.12	Animal Mode	32
3.12.1	Weigh Animal Using Manual Hold	32
3.13	Formulation Mode	33
3.13.1	Data Review	34
3.14	Unit Setting	34
3.15	Comparator Function	35
3.16	Addition Function	35
3.16.1	Weighing with Plus Side Addition	36
3.16.2	Weighing with Minus Side Addition	37
3.16.3	View/Delete Total Value	37
3.17	Tare-Subtraction Reminder Function	37
3.17.1	Set Tare-Subtraction Reminder Function	37
3.18	Zero-Point-Adjustment Reminder Function	38
3.18.1	Set Zero-Point-Adjustment Reminder Function	38
3.19	Stabilization Wait Function	38
3.19.1	Set Stabilization Wait Function	38
3.20	Bar Graph Display	39
3.21	Back Light Display	39
3.22	Auto Off Function	40
3.23	Simple Self Counting System (SCS) Method	40
3.23.1	Set SCS Method	40
4.0	Configuration	41
4.1	General Navigation	41
4.2	Numeric Value Entry	41
4.2.1	General Entry for All Models	41
4.2.2	For NTEP Only Models	41
4.2.3	For NTEP / Measurement Canada Models	41
4.3	Applications Menu	42
4.3.1	Operation Modes	43
4.3.2	NTEP Only Model Units	44
4.3.3	NTEP / Measurement Canada Model Units	45
4.3.4	Comparator Menu	46
4.3.5	Addition	47
4.4	Performance Menu	48
4.4.1	Stability	49
4.4.2	Response Speed	50
4.4.3	Zero Tracking	50
4.5	User Information	51
4.5.1	Preset Tare Mode	52
4.5.2	Input Preset Tare Value	52
4.5.3	Set the Discrimination Value of Comparator Function	53
4.6	External Input/Output Functions	54
4.7	Lock Functions	56
4.7.1	Total Lock Release	56
4.7.2	Key Lock Function	57
4.7.3	Menu Lock Function	57
4.8	Admin/Adjust Menu	58
4.8.1	Set Short Cut Mode	59
4.8.2	Set Free Keys	59



Technical training seminars are available through Rice Lake Weighing Systems. Course descriptions and dates can be viewed at www.ricelake.com/training or obtained by calling 715-234-9171 and asking for the training department.

4.8.3	Maintenance Settings	59
4.8.4	Balance Manage Menu	60
4.9	LCD Display Characters	61
4.9.1	7-Segment Characters	61
4.9.2	16-Segment Characters	61
5.0	Calibration	62
5.1	Perform External Calibration	62
5.2	Perform External Span Calibration Test	63
5.3	Sealing Instructions	64
5.3.1	Sealing Method 1	64
5.3.2	Sealing Method 2	65
6.0	Communications	67
6.1	RS-232 Connections	67
6.2	USB Connections	68
6.3	Basic Data Output Format	68
6.3.1	Data Description	69
6.4	CBM Data Output Format	70
6.4.1	Data Description	70
6.5	Input Commands	71
6.5.1	Input Command Composition 1	72
6.6	Command Formats	72
6.6.1	Input Command Composition 2	73
6.7	Response	74
6.8	External Contact Input	74
6.9	Communication Settings	75
6.10	Set Communication Parameters	75
6.10.1	Relay Contact Output (Option)	76
6.11	USB Communication and Bus Power Inputs	76
7.0	Troubleshooting and Maintenance	77
7.1	Maintenance Precautions	77
7.1.1	Cleaning—Round Pan Type	77
7.1.2	Cleaning—Square Pan Type	77
7.2	Error Messages	78
8.0	Specifications	81
8.1	Dimensions	81
8.2	Model Specifications	82
8.3	Basic Communication Specifications	82
8.4	Certifications and Approvals	82
8.5	Functional Specifications	83



Rice Lake continually offers web-based video training on a growing selection of product-related topics at no cost. Visit www.ricelake.com/webinars

1.0 Introduction

This document describes how to assemble and operate the Rice Lake TE Series Tuning Fork Enhanced Balance.



Manuals are available from Rice Lake Weighing Systems at www.ricelake.com/manuals

Warranty information is available at www.ricelake.com/warranties



Figure 1-1. Rice Lake TE Balances

Part No.	Model No.	Capacity	Readability (e)	Readability (d)	Pan Size	Weighing Chamber Size
187600	TE-223	220 g	0.01 g	0.001 g	4.6 in Diameter	8.9 x 5.9 x 5.9 in
186030	TE-623	620 g	0.01 g	0.001 g	4.6 in Diameter	8.9 x 5.9 x 5.9 in
186035	TE-3202	3,200 g	0.1 g	0.01 g	6.3 x 7.1 in	—
186036	TE-6202	6,200 g	0.1 g	0.01 g	6.3 x 7.1 in	—
186037	TE-15001	15,000 g	1 g	0.1 g	6.3 x 7.1 in	—

Table 1-1. Available NTEP Only Models

Part No.	Model No.	Capacity	Readability (e=d)	Pan Size	Weighing Chamber Size
204688	TE-322NC	320 g	0.01 g	4.6 in Diameter	8.9 x 5.9 x 5.9 in
204689	TE-1501NC	1500 g	0.1 g	6.3 x 7.1 in	—
204690	TE-8200NC	8200 g	1 g	6.3 x 7.1 in	—

Table 1-2. Available NTEP/Measurement Canada Models

Part No.	Description
186074	Power supply 100-240 VAC, 50-60 Hz
186075	In-use dust cover 1,200 g models and lower
186076	In-use dust cover 3,200 g models and higher
186077	Specific Density Measurement Kit
186079	Ethernet TCP/IP option card

Table 1-3. Available Options

1.1 Safety

Safety Definitions:



DANGER: Indicates an imminently hazardous situation that, if not avoided, will result in death or serious injury. Includes hazards that are exposed when guards are removed.



WARNING: Indicates a potentially hazardous situation that, if not avoided, could result in serious injury or death. Includes hazards that are exposed when guards are removed.

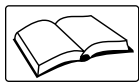


CAUTION: Indicates a potentially hazardous situation that, if not avoided, could result in minor or moderate injury.



IMPORTANT: Indicates information about procedures that, if not observed, could result in damage to equipment or corruption to and loss of data.

General Safety



Do not operate or work on this equipment unless this manual has been read and all instructions are understood. Failure to follow the instructions or heed the warnings could result in injury or death. Contact any Rice Lake Weighing Systems dealer for replacement manuals.



WARNING

Failure to heed could result in serious injury or death.

Do not allow minors (children) or inexperienced persons to operate this unit.

Do not operate unless unit is completely assembled.

Do not use for purposes other than weight taking.

Do not exceed the rated specification of the unit.

Do not remove or obscure warning labels.

Do not use in wet locations or with wet hands.

Do not use solvents or aggressive substances to clean the unit.

Review MSDS (Material Safety Data Sheets) when applicable.

Only use the specified power supply supplied by Rice Lake Weighing Systems.

Do not disassemble or modify the unit. For inspection and adjustments, contact Rice Lake Weighing Systems.

Do not use in an explosive environment.

1.2 FCC Compliance

United States

This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is likely to cause harmful interference in which case the user will be required to correct the interference at his own expense.

Canada

This digital apparatus does not exceed the Class A limits for radio noise emissions from digital apparatus set out in the Radio Interference Regulations of the Canadian Department of Communications.

Le présent appareil numérique n'émet pas de bruits radioélectriques dépassant les limites applicables aux appareils numériques de la Class A prescrites dans le Règlement sur le brouillage radioélectrique édicté par le ministère des Communications du Canada.

1.3 Ensuring Accurate Measurements

To maintain consistent and accurate measurements, eliminate error-factors by following the precautions below.

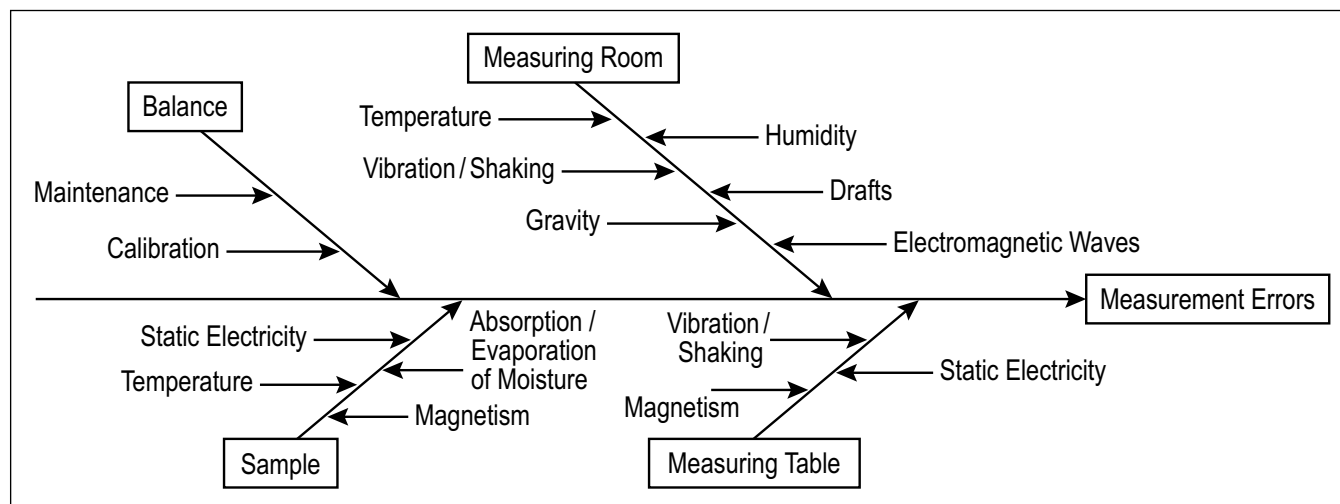


Figure 1-2. Measurement Error Factors

1.3.1 Measuring Environment Precautions

Temperature / Humidity / Air Draft / Atmospheric Pressure / Dust

- Keep the room temperature constant to avoid condensation and indication drift
- Avoid exposure to direct sunlight
- Avoid exposure to airflow (air conditioner, heat ducts) and atmospheric pressure on the balance
- Avoid measurement in dusty locations and areas with low humidity that can cause static electricity

Vibration / Shaking

- When possible, locate balance in a room on the first floor or basement
- Rooms near a road/railroad should be avoided
- Place balance on a table/counter that is not affected by vibration
- Do not place cloth or paper under the balance

Gravity

- The latitude and altitude of a measuring location can affect a weight reading due to changes in gravity
- Calibrate the balance at a measuring location to account for location gravity

Electromagnetic Wave

- Avoid locations where strong electromagnetic wave generating objects are present
- Avoid using tables that are subject to magnetism or static electricity

1.3.2 Specimen Related Precautions

Static Electricity

- Weighing electrically charged specimens make the displayed value unstable, reducing the reproducibility of a test result
- Neutralize electrically charged specimens before measurement, such as synthetic resin or glass made materials

Magnetism

- Eliminate the magnetism or place the specimens on a container/plate to protect the weighing mechanism

Moisture Absorption/Evaporation

- Put wet or volatile specimens in a container with a tight fitting lid to prevent evaporation during the weighing process

Specimen Temperature

- If specimen temperature is excessively high or low, allow its temperature to stabilize at room temperature prior to measuring
- Ensure the draft shield interior temperature is equal to room temperature prior to measuring
- Body temperature of operator affects measurement results; handle specimens with tweezers and refrain from putting hands in the draft shield during operation

1.3.3 Balance Main Unit Precautions

Operating precautions

- If a dust cover is used, wipe with an anti-static agent or remove the cover
- For stability, turn on the balance 30 minutes or more and load the balance a few times with a calibration weight equivalent to the weighing capacity prior to using

Adjustment

- Calibrate periodically with a certified test weight with a finer resolution than the balance
- Adjustment is needed when:
 - Using the balance for the first time
 - Using the balance after a long period of non-use
 - Relocating the balance
 - A major change in temperature, humidity or atmospheric pressure has happened

Maintenance

- Remove any dust or liquid from the pan and/or pan base prior to operation
- Ensure that no dust or liquid enters the balance when cleaning
- Frequent cleaning of the balance is required (see [Section 7.1 on page 77](#))

2.0 Installation



IMPORTANT: Use balance only in areas free from environmental conditions that could affect the accuracy ([Section 1.3 on page 9](#)). Legal for Trade balance must be verified and sealed ([Section 5.3 on page 64](#)).

2.1 Unpack Balance

Unpack the balance and inspect the contents. Report missing or damaged components to the shipper and Rice Lake Weighing Systems immediately.

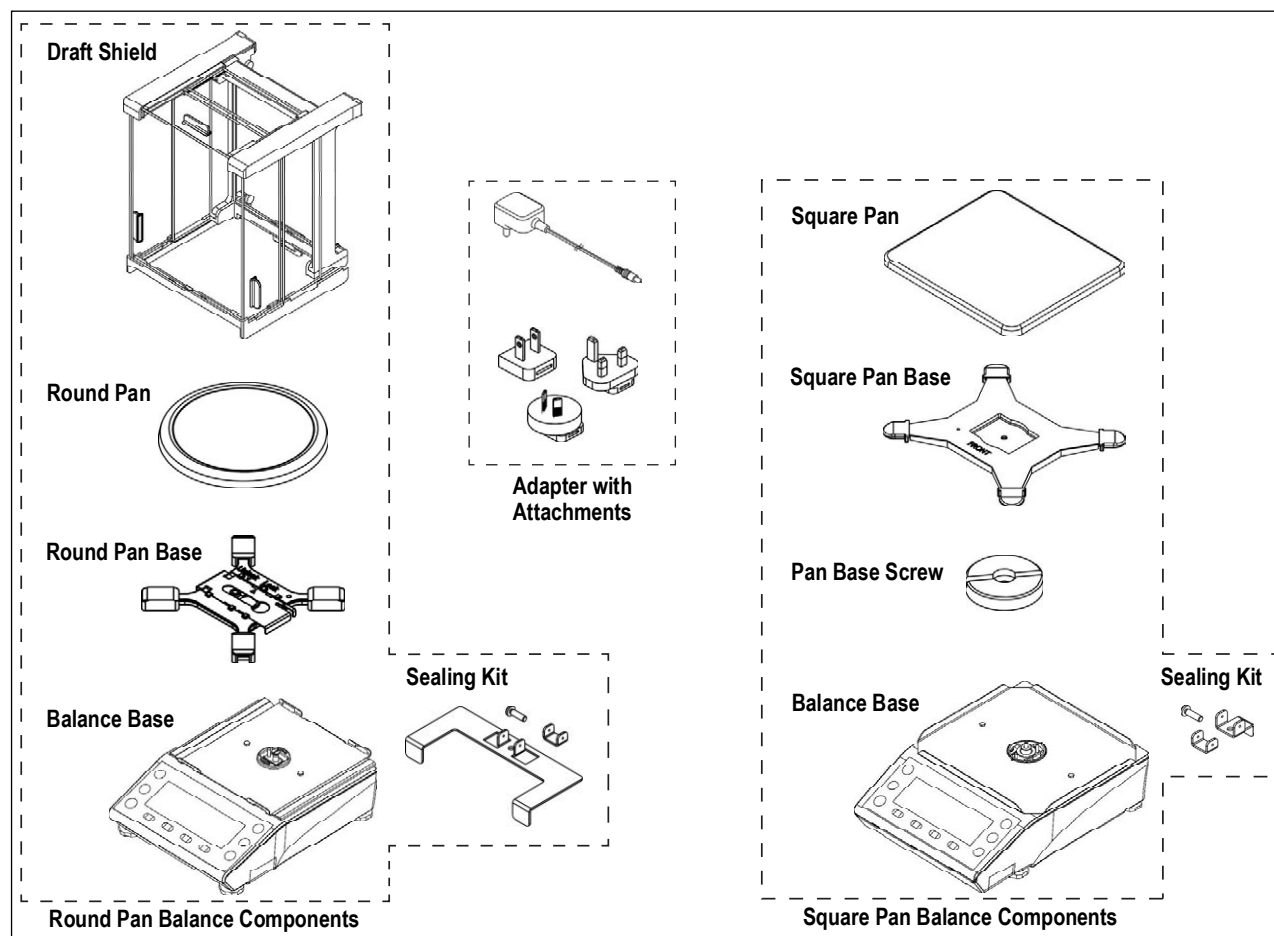


Figure 2-1. Balance Components

2.1.1 Connections and Part Locations

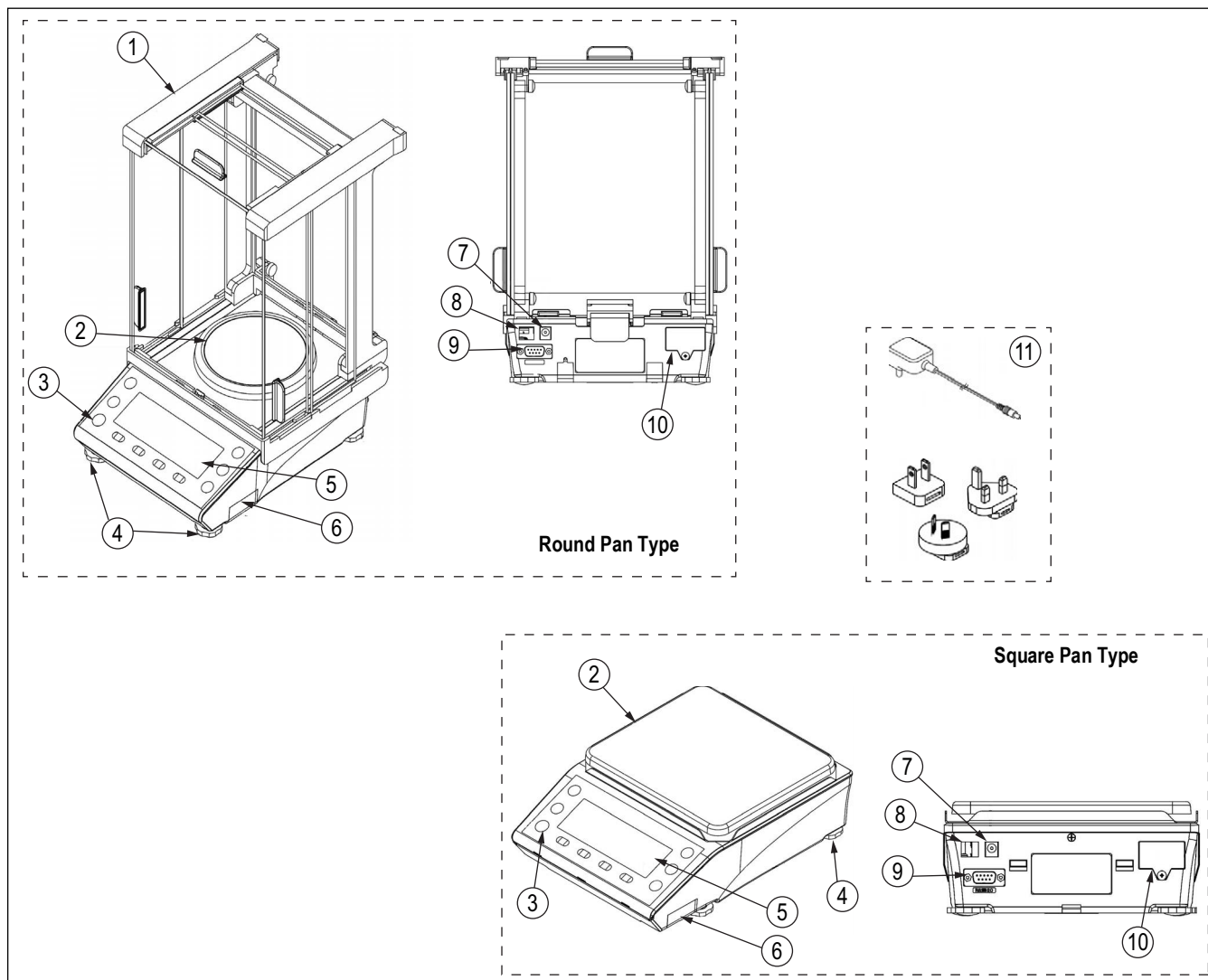


Figure 2-2. Connection and Part Locations

Item No.	Description
1	Draft shield (Section 2.5 on page 16)
2	Weighing pan
3	Level
4	Adjustable Feet
5	Display
6	Battery case
7	AC adapter jack
8	USB connector (Type B)
9	RS-232C connector (D-sub 9 pin male)
10	Option slot
11	Adapter with attachments

Table 2-1. Parts and Locations

2.2 Assemble Balance

Use the following steps to assemble the balance.

2.2.1 Round Pan Base Assembly

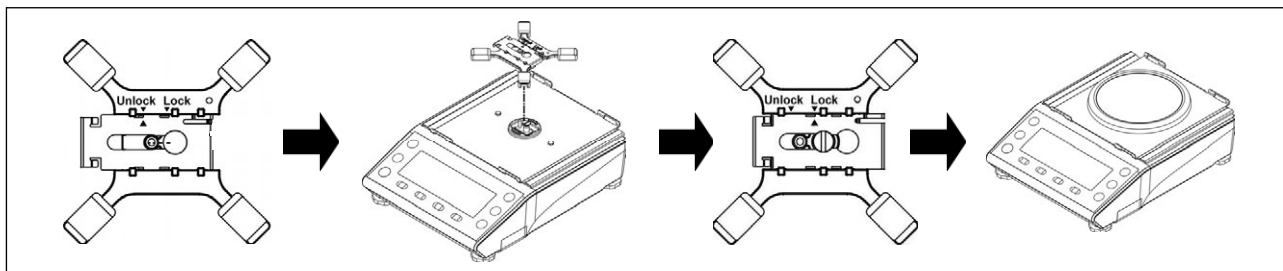


Figure 2-3. Assemble Base

1. Ensure the lock plate is in the unlocked position.
2. Set the pan base on the balance.
3. Slide the lock into the locked position.
4. Mount the weighing pan.
5. Connect the AC adapter or insert batteries ([Section 2.3 on page 14](#)).
6. Assemble and install the draft shield ([Section 2.5 on page 16](#)).

2.2.2 Square Pan Base Assembly

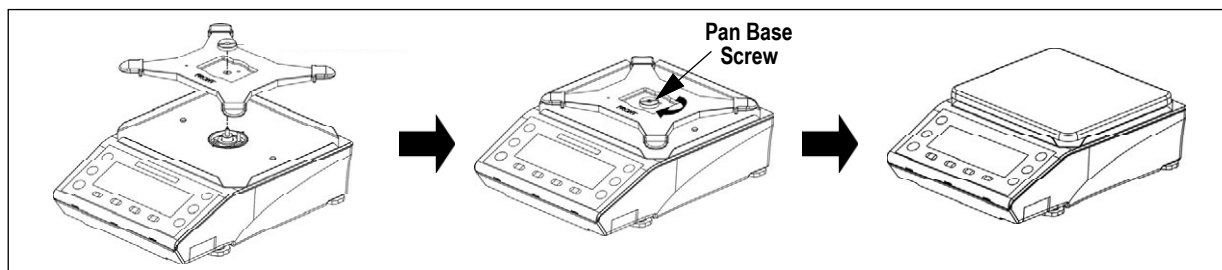


Figure 2-4. Square Base Assembly

1. Ensure the word **Front** is aligned towards the display of the balance and attach the pan base to the balance.
2. Tighten the pan base screw firmly.
3. Place the weighing pan on the pan base.
4. Connect the AC adapter or insert batteries ([Section 2.3 on page 14](#)).

2.3 Battery Installation

Use the following steps to install four AA batteries into the unit. Alkaline, manganese and nickel-metal hydride batteries can be used. Approximate battery life is 150 hours using alkaline batteries with the back-light and external output off.

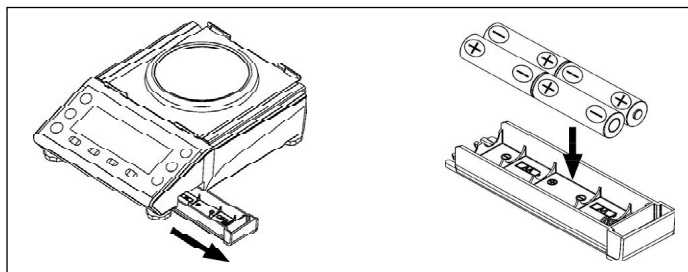


Figure 2-5. Install Batteries

1. Pull out the battery case.
2. Insert four AA batteries into battery case. Make sure to insert batteries with the positive and negative poles correctly inserted.
3. Insert battery case back into the unit until it clicks in place. When the balance is battery operated, one of the icons below will display:

 Full Charge

 Low Charge

 Replace Batteries



IMPORTANT: Observe the following statements when changing or replacing the batteries.

- If the balance is not going to be used for a long period of time remove the batteries
- Dispose of used batteries in accordance with state and local regulations
- Do not use batteries that leak or are damaged
- Only use new batteries in the balance, never mix used and new batteries or different brands/manufacturers

2.4 Leveling the Balance

Use the following steps to level the balance.

1. Turn the adjustable feet as shown in [Figure 2-6](#) to unlock them.

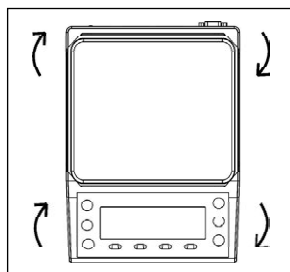


Figure 2-6. Unlock Adjustable Feet

2. Adjust balance feet until the bubble level is centered in the center circle. See [Figure 2-7](#) for illustrated directions for adjusting the feet to level the balance.

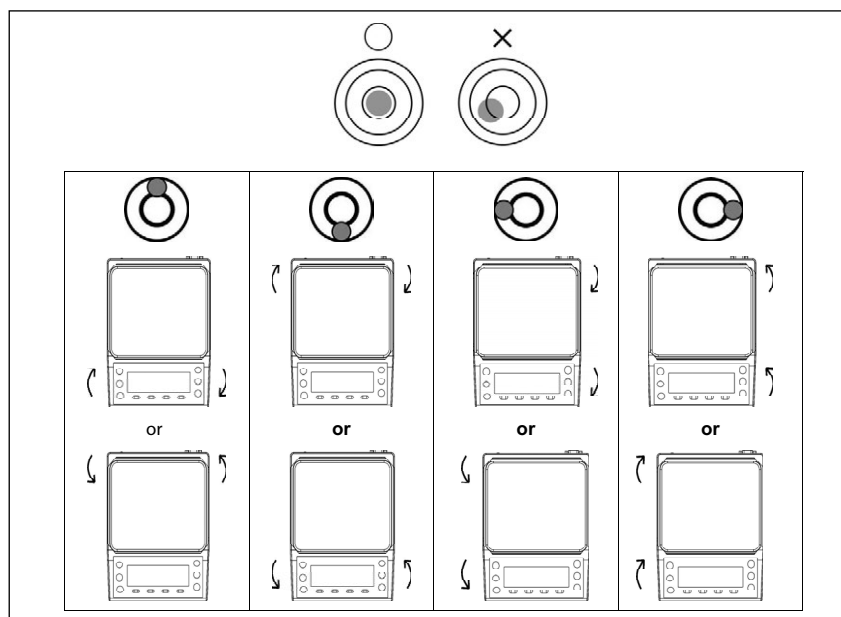


Figure 2-7. Level the Balance

2.5 Draft Shield Assembly

The draft shield for round pan balance models must be assembled using the following steps.

1. Place the draft shield base on a flat, level surface.
2. Insert one of the rear panel supports into the bottom frame as shown in [Figure 2-8](#). Ensure flat side is on the outside.
3. Secure the rear panel support to the base with a provided screw and repeat for the other side.

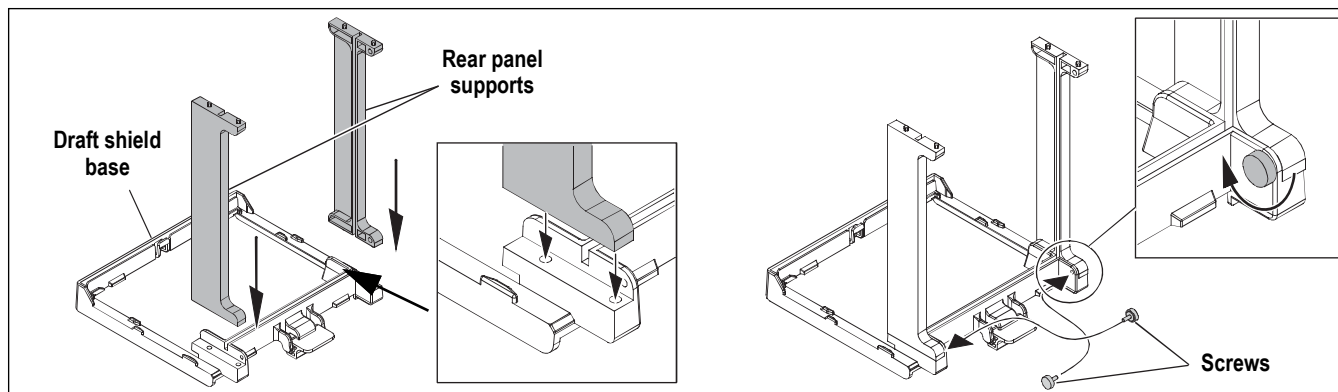


Figure 2-8. Install Rear Panel Supports

4. Insert the rear glass panel into the slots of the back supports.
5. Gently place top assembly support on the rear panel supports.

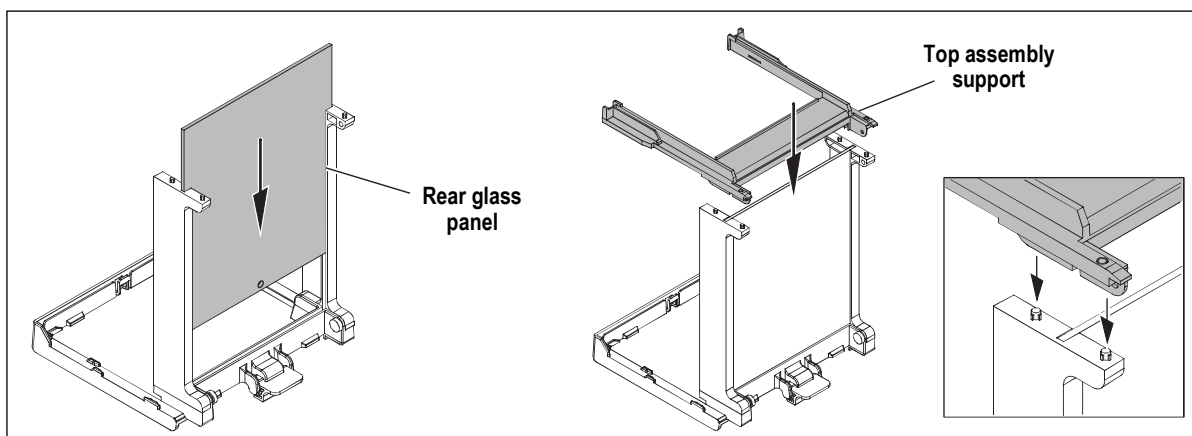


Figure 2-9. Insert Front and Rear Panels

6. Insert the front glass panel into the base, ensuring it clicks as shown in [Figure 2-10](#).
7. Carefully connect to the front glass panel.

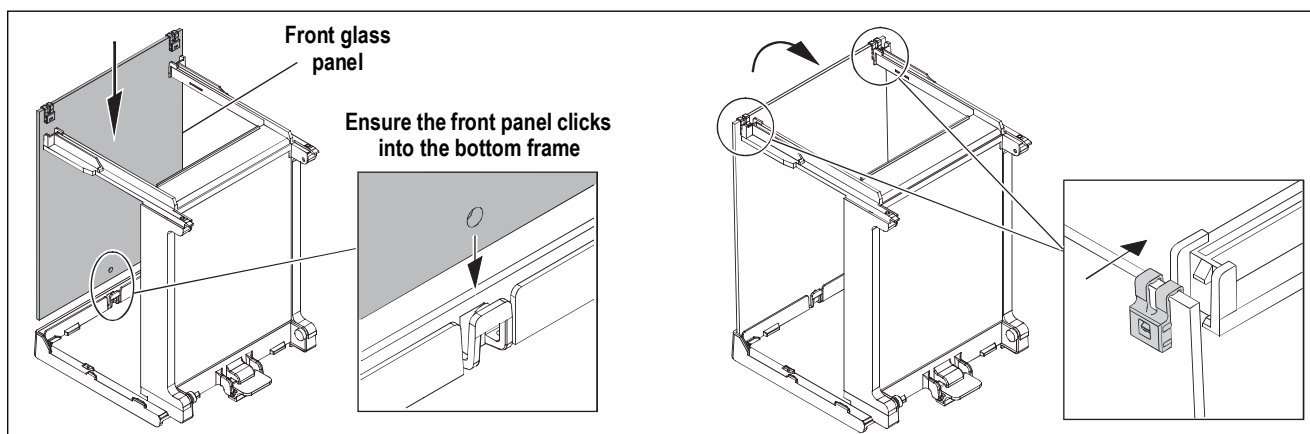


Figure 2-10. Install Top Assembly

8. Place the upper side brackets on top of the assembly as shown in [Figure 2-11](#) and secure with provided screws.

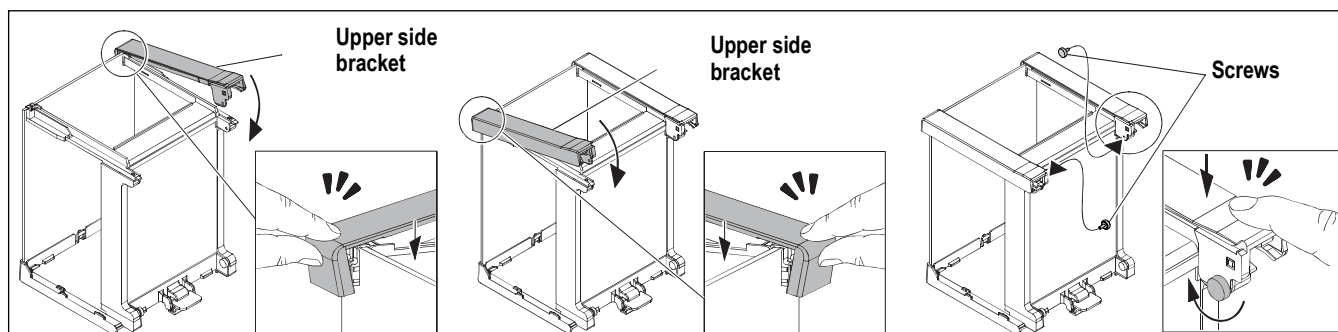


Figure 2-11. Install Upper Brackets

9. Slide the top door assembly, with the handle towards the front, into the upper brackets as shown in [Figure 2-12](#). Ensure the side rails of the top assembly are in the grooves in the top brackets.

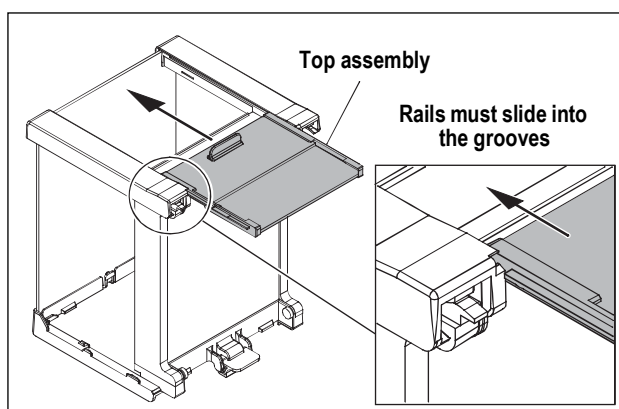


Figure 2-12. Install Top Panel

10. Place the glass side panels together ensuring the rollers are on the outside. Use one panel with the handle and one without as a pair.
11. Hold the panels together, ensuring the panel with the handle is facing outward. With the panels held tightly together, insert the rollers into the frame of the top assembly as shown in [Figure 2-13](#).

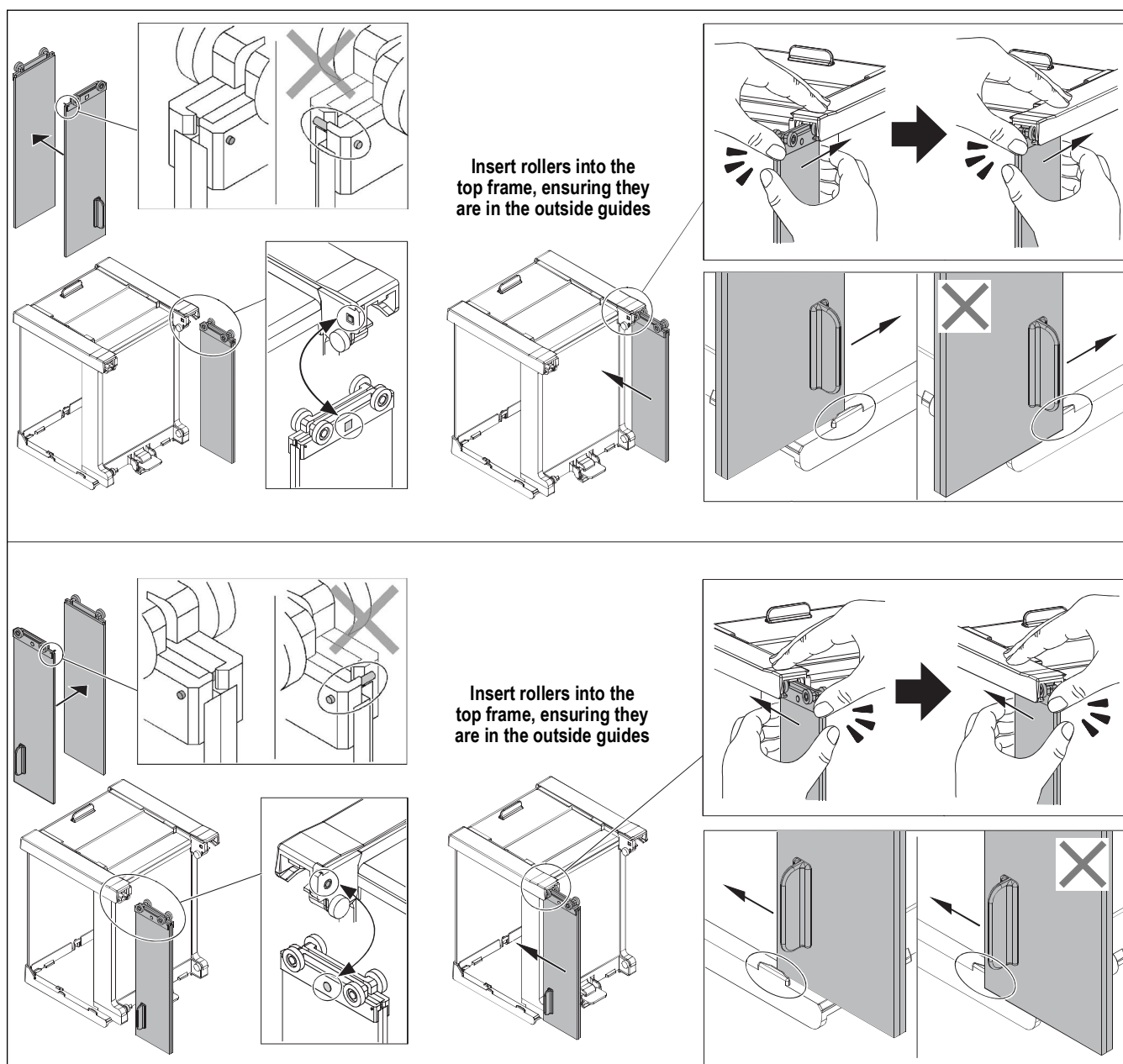


Figure 2-13. Insert Side Assemblies

12. Push the pair as far as they will go, there is a stop for the inside panel.
13. Slide the outer panel to the front of the assembly.

14. Ensure the latch on the back of the draft shield is in the up position.
15. Place the completed draft shield assembly on the balance as shown in [Figure 2-14](#).

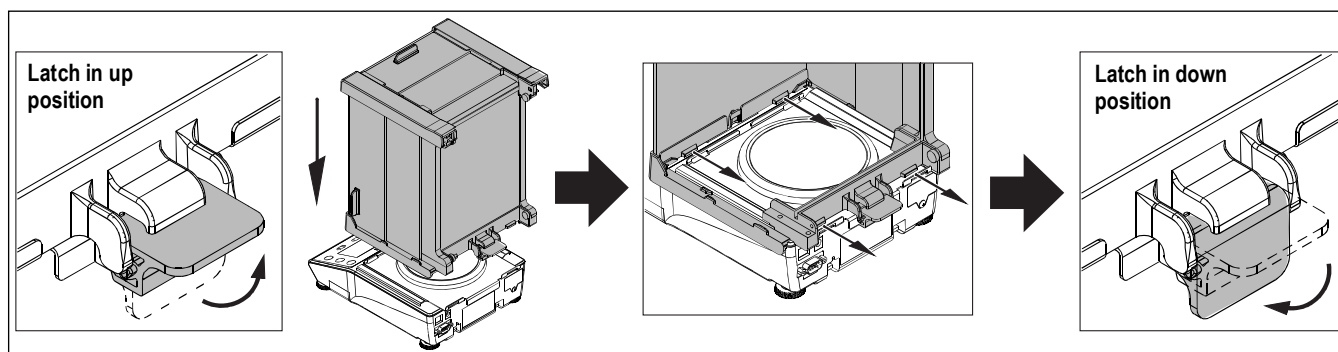


Figure 2-14. Install Draft Shield on Balance

16. Gently position the draft shield assembly toward the back of the balance to seat it in the slots of the balance.
17. Rotate the latch to the down position securing the draft shield assembly to the balance.

! IMPORTANT: Do not move balance by holding the draft shield. Make sure to hold the main body of the balance to move it.

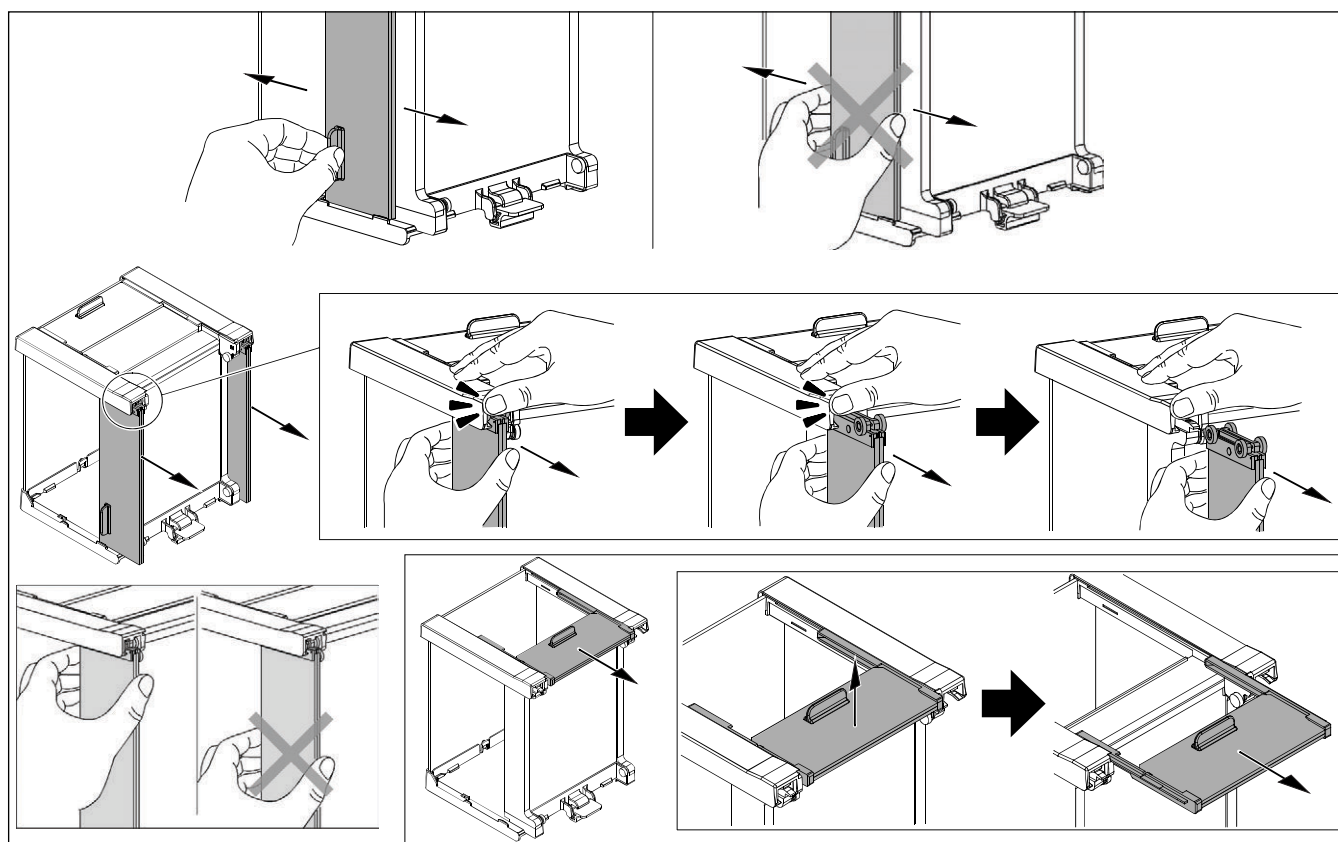


Figure 2-15. Use and Disassembly

2.5.1 Draft Shield Assembly Parts

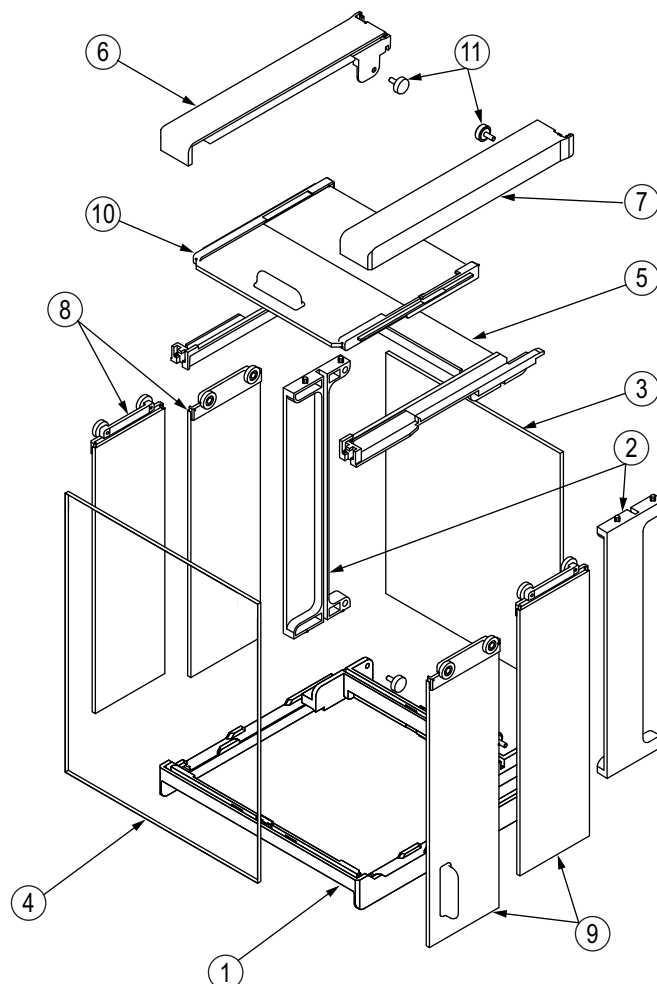


Figure 2-16. Draft Shield Parts

Item No.	Description
1	Draft Shield Base
2	Rear Panel Support, Plastic
3	Rear Glass Panel
4	Front Glass Panel
5	Top Assembly Support
6	Left Upper Side Bracket
7	Right Upper Side Bracket
8	Left Door Assembly
9	Right Door Assembly
10	Top Door Assembly
11	Screw (No. 1, White, d=3, L=8)

Table 2-2. Draft Shield Parts List

3.1.3 Front Panel Item Descriptions

Item No.	Description
1	Power: Press to turn unit on; a long press to turn off
2	Menu: Press to enter/exit the setting menu; cancels setting value selection; returns to the operation mode
3	Output: Used for data outputting; data importing in the Statistics/Formulation mode
4	Zero: Zeros the balance
5	Tare: Allows the container/vessel weight to be removed from the gross weight, leaving the net weight displayed; only available for NTEP only models
6	F1: Use for selecting mode, function or item; use for moving up to menu/item selections; use for incrementing numeric values
7	F2: Use for selecting mode, function or item; use for moving down to menu/item selections; use for decrementing numeric value
8	F3: Select mode, function or item; use for moving to the upper menu layer; or use for selecting the digit to change
9	F4: Select mode, function or item; enter a menu; select digit to change; returning to the setting menu/weighing mode
10	Bubble Level: Indicates when the balance is level and ready for use
11	Asterisk: Lights in the standby status; indicates addition available status when the adding function is used
12	Negative: Indicates the negative weight value and numeric
13	Stable: Indicates balance is stable when displayed
14	Zero: Indicates the balance is at zero
15	Battery Level: Indicates level of battery power available
16	Colon: used when date and time display
17	Bar Graph: Indicates the present total amount relative to the weighing capacity defined as 100%; indicates the state of span adjustment/calibration with internal weight
18	7-Segment Main String: indicates the weight value; indicates the simplified character
19	Auxiliary Balance Interval: Lit when the auxiliary balance interval is displayed; not available in Legal for Trade models
20	Animal Weigh Mode: Displays when in Animal Weigh Mode; not available in Legal for Trade models
21	Output: Displayed when data is being output to external devices
22	Brutto/Gross: Indicates gross weight is being displayed Net: Indicates a tare weight is being subtracted; indicates a preset tare weight; Only available for NTEP only models Preset Tare: Indicates a preset tare weight; Only available for NTEP only models
23	16-Segment Unit String: Indicates various units
24	Units: Indicates unit of value displayed; mg is not available in Legal for Trade models
25	16-Segment Message String: Displays various messages


Table 3-1. Front Panel Item Descriptions




NOTE: The functions of the F-keys ↑, ↓, ←, →, ↵, or ▼ appear on the display directly above the corresponding F-key. Shortcuts for various modes/functions can be assigned to F-keys (Section 4.8 on page 58).

3.2 Power On/Off

When balance is plugged in, ✱ displays.



Press . A self-check runs, then 000000 displays.

- Do not press any buttons during the self-check
- The balance starts in the last operation mode set; Legal for Trade balances always open in Weigh Mode

Press and hold  to put the balance into standby, ✱ displays.

3.2.1 Standby

When in standby, ***** displays.





- Press . A self-check runs and balance enters operation mode
- Do not press any buttons during the self-check
- Press and hold  to place the balance into standby during operation
- The balance starts in the operation mode it was in before it went into standby
- If using battery power, ***** is not displayed in standby

3.2.2 Balance Operation Check

Press on the weighing pan lightly to see if the weight display changes.

3.3 Display Navigation

Use the F-Keys to navigate through the menus and enter numerical values.

-  = ↑ Increments numeric values. Scrolls through menus/selections
-  = ↓ Decrements numeric value. Scrolls through menus/selections
-  = ← Select digits and returns to previous level
-  = → or ↵ Enters the value or a level



NOTE: The functions of the F-keys ↑, ↓, ←, →, ↵, or ▼ appear on the display directly above the corresponding F-key.

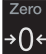
3.4 Zero-Point Adjustment

A zero-point adjustment cannot be performed if the weight is over the zero-point adjustment range.




NOTE: On Legal for Trade balances, this setting is read only, the wait time cannot be changed.

Use the following steps to do a zero-point adjustment.

1. Ensure the weigh pan is empty.
2. Press , 000000 displays.

The stability wait time can be set in the Applications menu, see [Section 4.3 on page 42](#).

3.4.1 Zero-Point Adjustment Limits

The **Zero Point** limit is the threshold at which the operator can press  and have it zero. There is an upper and lower threshold. The balance cannot be zeroed outside of that range. The zero range is 1.5% of capacity above or below zero.

Model	Lower Limit (g)	Upper Limit (g)
TE-223	-3.300	3.300
TE-623	-9.300	9.300
TE-3202	-48.00	48.00
TE-6202	-93.00	93.00
TE-15001	-225.00	225.00

Table 3-2. Zero-Point Adjustment Limits (NTEP Only Models)

Model	Lower Limit (g)	Upper Limit (g)
TE-322NC	-4.80	320
TE-1501NC	-22.50	1500
TE-8200NC	-123.00	8200

Table 3-3. Zero-Point Adjustment Limits
(NTEP/M Measurement Canada Models)

3.5 Tare Value



NOTE: The tare function is only available for NTEP only models.
NTEP/Measurement Canada models do not support the tare function.

To weigh using a container, the container weight should be subtracted from the total value. Then results displayed are the weight of the product only.



NOTE: A Tare Value is included in the total maximum capacity.
When turning on the power, placing a tare that exceeds the zero point adjustment range at the time of power supply, the tare subtraction is executed.

1. Place the empty container on the balance. The weight of the container displays.
2. Press . 000000 and the **Net** icon displays.
3. Place the product to be weighed in the container. The net weight displays.
4. Remove the product and container from the balance.
5. Press to remove the tare.

3.5.1 Preset Tare

When a tare weight is already known, the tare subtraction can be configured as a **Preset Tare** parameter. Five preset tare values can be stored. See [Section 4.5.2 on page 52](#) for steps to assign values to the tare presets 1-5.

To select a preset tare:

1. Press . displays.
2. Press or to scroll to **USER INFO**.
3. Press . **31 PT MODE** and the current setting displays.
4. Press . The current setting begins to flash.
5. Press or to select desired preset tare number (1-5).
6. Press to save the setting.
7. Press to return to operation display. **NET PT** displays in the upper right.

To exit a preset tare:

1. Ensure there is no weight on the balance.
2. Press . **ZERO** flashes and display returns to operation mode. **NET PT** is no longer displayed.

3.5.2 Check Tare Weight Using an F-Key

If an F-Key has been set to Tare, it can be used to see what the tare weight is.

1. Set an F-Key to Tare. See [Section 4.8 on page 58](#).
2. With tare in place, press until **TARE** is displayed above an F-Key.
3. Press the **TARE** F-Key. The tare weight displays.
4. Press to return to operation mode.

3.5.3 Add to Product

To weigh additional product without the value of the existing product, use the following steps.

1. Place first product to be weighed on the balance.
2. Push . 000000 displays.
3. Add additional product to the balance/container. The weight of additional product only is displayed.

3.6 Weigh Mode

Weigh mode is the basic mode for weighing product. Place the product on the tray, the weight will display in the units set in configuration. Each mode of operation has F-Key commands available.



NOTE: Weigh mode can be customized for what is displayed. The F1-F3 keys on the primary and secondary pages of the display can be any selection from the 62 Free Key menu ([Section 4.8 on page 58](#)).

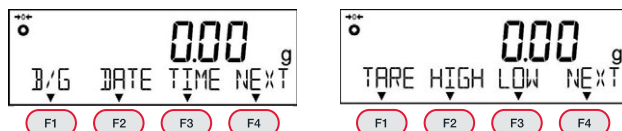


Figure 3-3. Weigh Mode F-Key Command Defaults (NTEP Only Models)

Command	Description
B/G	Toggles Brutto/Gross command ON and OFF, showing gross when ON or net (if a tare exists) when OFF; only available for NTEP only models
DATE/TIME	Date and Time commands display current date or time respectively when ON or pressed
NEXT	Press to move between menu or F-key selections
TARE	Displays current tare; only available for NTEP only models
HIGH/LOW	High and Low commands display current high or low limit and allows new entry of high or low limit respectively
g/mg/ct	Unit of measure commands switch weight display to unit set (g = grams / mg = milligrams / ct = carats)
ICAL/CAL	Calibration commands display current internal span adjustment (ICAL) or external span adjustment (CAL); ICAL only available for NTEP only models
HOLD	Toggles Hold command ON and OFF, when ON the weight measurement indication holds on the display

Table 3-4. Weigh Mode F-Key Commands



NOTE: Refer to the 61 Short Cut Mode menu ([Section 4.8 on page 58](#)) to set F1-F3 keys for mode selection options.

3.7 Counting Mode

Counting Mode counts the number of items placed on the balance. There are two methods to input the unit weight.

- Actual Value Setting Method: Place the specified number of samples on the balance to record the average unit weight
- Numeric Value Setting Method: Input numeric value of the unit weight by key operation

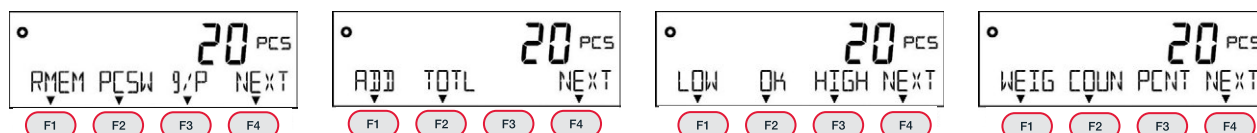


Figure 3-4. Counting Mode F-Key Commands

Command	Description
RMEM	Change the unit weight
PCSW	Unit weight
g/P	Weight of samples/number of samples flash on display
ADD	Execute addition
TOTL	Total value
LOW	Displays onW or NUM. Select onW to display current limit low, press NUM to enter a new limit low
OK	Displays onW or NUM. Select onW to display current limit OK, press NUM to enter a new limit OK
HIGH	Displays onW or NUM. Select onW to display current limit high, press NUM to enter a new limit high
WEIG	Press to display Weigh Mode
COUN	Press to display Count Mode
PCNT	Press to display Percentage Mode
NEXT	Press to move between menu selections


Table 3-5. Counting Mode F-Key Commands

3.7.1 Actual Value Setting

Place the specified number of product on the balance to record the average unit weight internally.

Press **F3** or **F4** to select No (change) or Yes (don't change). This step is not available if no data is recorded.

For YES (No Change):


1. Press **F4** (**YES**).
2. Press **F4** to scroll through the F-Key commands until **OK** is displayed.
3. Press the **YES** F-Key.
4. Place a container on the weigh pan and push  if needed. The zero-point adjustment or tare is set.



NOTE: For NTEP/Measurement Canada models, press  to set the zero-point adjustment.

5. Add product for counting.

For NO (Change):

1. Press **F3** (**NO**). Press **F1** or **F2** to scroll through selections.
 - ON 5: 5 PCS
 - ON 10: 10 PCS
 - ON 30: 30 PCS
 - ON 50: 50 PCS
 - ON 100: 100 PCS
 - ON VAR: 1-999 PCS
 - PCSWGT: Unit weight value input
2. With the desired selection displayed, press **F4**.
3. Place a container on the weigh pan and push  if needed. The zero-point adjustment or tare is set.




NOTE: For NTEP/Measurement Canada models, press  to set the zero-point adjustment.

4. Add product for total count.

3.7.2 Numeric Value Setting Method

Use the key operation to input a numeric value.

1. Press **F3** or **F4** to select whether or not to use the previous data.
When there is no data record, this step is skipped.
2. Press **F1** or **F2** to scroll through F-Key commands until **YES** displays. When **YES** is selected, skip to Step 5.
3. Select the unit weight value input mode by pressing **F1** or **F2** to PCSWGT.
4. Use **F1** or **F2** to enter the weight/number of the product.
5. Press **F4** to fix and the unit weight is recorded.
6. Place a container on the weighing pan and press .



NOTE: Press  for NTEP/Measurement Canada models.

7. Place samples on the weighing pan and the count result is displayed.

3.8 Percentage Mode

The weight of a product to be weighed is shown in a percent relative to the reference weight.

There are two methods to enter the reference weight:

- Actual Value Setting Method [onW]: Where placing the reference weight on the balance is done to record the weight.
- Numeric Value Setting Method [NUM]: Inputting the numeric value of the reference weight is done by a key operation.

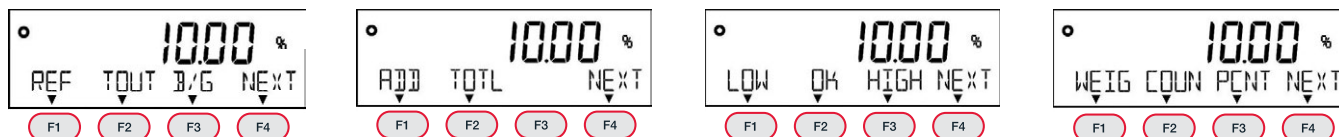


Figure 3-5. Percentage Mode F-Keys (NTEP Only Models)



Command	Description
REF	Display reference value
TOUT	Tare value output; only available for NTEP only models
B/G	Gross; only available for NTEP only models
ADD	Executes addition
TOTAL	Total value
LOW	Displays onW or NUM. Select onW to display current limit low, press NUM to enter a new limit low
OK	Displays onW or NUM. Select onW to display current limit OK, press NUM to enter a new limit OK
HIGH	Displays onW or NUM. Select onW to display current limit high, press NUM to enter a new limit high
WEIG	Press to display Weigh Mode
COUN	Press to display Count Mode
PCNT	Press to display Percentage Mode
NEXT	Press to move between menu selections

Table 3-6. Percentage Mode F-Key Commands



NOTE: ADD and TOTL can be used when 14 ADDITION is set to Valid. See [Section 4.3 on page 42](#).

Adding function display is not available on verified Legal for Trade balances.

1. Press , then press **F4** to view the current operation mode.
2. If needed, press **F4** and use **F1** or **F2** to scroll to **PCNT**.
3. Press **F4**, then press  to return to operation mode.
4. To use the previous data, press **F4** (**YES**). Skip to Step 6. (When there is no data record, this step is skipped) or
To set the percent, press **F3** (**NO**).
5. Press **F3** or **F4** to set one of the following:
 - Actual value (**onW**): Placing the reference weight on the balance and press **F4** (**OK**)
 - Numeric value (**NUM**): Use **F1** or **F2** to enter the reference weight and press **F4** (**OK**)
6. Weigh the product. The ratio of the weight of the product to the reference weight is displayed in percent.

3.9 Coefficient Mode

Measured weight is multiplied by the preset coefficient. This mode is not available in Legal for Trade units.

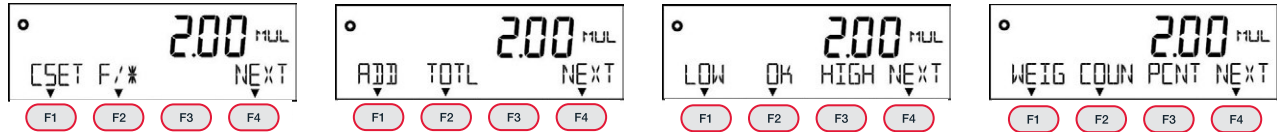




Figure 3-6. Coefficient Mode F-Keys

Command	Description
CSET	Set the Coefficient value <ul style="list-style-type: none"> NO: Allows the change of the current coefficient YES: Accept the currently set coefficient
F/*	Toggles between MUL and g
ADD	Executes addition
TOTAL	Displays total value
LOW	Displays onW or NUM. Select onW to display current limit low, press NUM to enter a new limit low
OK	Displays onW or NUM. Select onW to display current limit OK, press NUM to enter a new limit OK
HIGH	Displays onW or NUM. Select onW to display current limit high, press NUM to enter a new limit high
WEIG	Press to display Weigh Mode
COUN	Press to display Count Mode
PCNT	Press to display Percentage Mode
NEXT	Press to move between menu selections

Table 3-7. Coefficient Mode F-Key Commands



NOTE: ADD and TOTL can be used when 14 ADDITION is set to valid. See [Section 4.3 on page 42](#).

1. Press , then press **F4**. The current operation mode displays.
2. If needed, press **F4**. The currently displayed mode flashes.
3. Use **F1** or **F2** to scroll to **MULT**.
4. Press **F4**, then press  to return to operation mode.
5. To use the existing coefficient value, press **F4** (**YES**). (When there is no data record, this step is skipped.)
or
Press **F3** (**NO**) to set the coefficient value. See [Section 3.3 on page 23](#).
6. Weigh the product. The weight is multiplied by the coefficient and the result displays.
The coefficient can be changed at anytime using the **CSET** F-Key.

3.10 Specific Gravity Mode

In the Specific Gravity Mode, the ratio of the density of a substance to the density of water at its densest (4° C for liquids) is calculated.

3.10.1 Menus Available in Operation Mode

Use F-Keys to select and scroll through menu items.



Figure 3-7. Specific Gravity Mode F-Keys

Command	Descriptions
RSET	Select the liquid <ul style="list-style-type: none"> • OTHER: Any liquid other than water • H2O: Using water
WAIR	Weight in air displays
WLIQ	Weight in liquid displays
WEIG	Press to display Weigh Mode
COUN	Press to display Count Mode
PCNT	Press to display Percentage Mode
NEXT	Press to move between menu selections or measuring mode switching

Table 3-8. Specific Gravity Mode F-Key Commands

3.10.2 Materials Required in Specific Gravity Mode



Materials needed to measure a specific gravity include:

- Water tank
- Hanging string or wire
- Basket for placing the sample
- Thermometer

3.10.3 Prepare Equipment to Measure Specific Gravity

1. Input the water temperature or the specific gravity of the reference liquid.
2. Measure the sample weight in the air.
3. Compensate the buoyancy acting on the basket.
4. Measure the sample weight in the liquid.
5. The specific gravity of the sample is displayed.

3.10.4 Measure the Specific Gravity Mode

1. Press , then press **F4**. The current operation mode displays.
2. If needed, press **F4**. The currently displayed mode flashes.
3. Use **F1** or **F2** to scroll to **SPGR**.
4. Press **F4**, then press  to return to operation mode.

5. Select the reference liquid by pressing **F3** or **F4**;
 - OTHER: Liquid other than water
 - H2O: Water
6. Enter the specific gravity of the reference liquid and press **F4** to save.
7. Set the net/basket on the balance and press **Tare**.

 **NOTE:** Press **Zero** for NTEP/Measurement Canada models.

8. Add material/liquid to the net/basket to measure the weight.
9. Press **F4** to save.
10. Remove the material/liquid on the net/basket and press **Tare**.

 **NOTE:** Press **Zero** for NTEP/Measurement Canada models.

11. Lower the net/basket into the liquid.
12. Press **Tare** to compensate the buoyancy acting on the net/basket.

 **NOTE:** Press **Zero** for NTEP/Measurement Canada models.

13. Put the material/liquid on the basket into the liquid and press **F4** to save. The specific gravity of the specimen is automatically calculated and displayed.

3.11 Statistics Mode

The statistical operation function collects weight data and indicates maximum, average and other statistical values. Not available in Legal for Trade balances.

- Only mg or g can be used
- Each calculation result, except CV, follows the smallest readability used to record the weighing data
- Up to 999 weight data can be saved

3.11.1 Menus Available in Operation Mode

Use F-Keys to select and scroll through menu items.






Figure 3-8. Statistics Mode F-Keys

Command	Description
POUT	Print
DATA	Display data menu options
DEL	Delete data <ul style="list-style-type: none"> • ALL: Deletes all stored data for Statistics Mode • LAST: Deletes only last data saved • RET: Return to operation mode
WEIG	Press to display Weigh Mode
COUN	Press to display Count Mode
PCNT	Press to display Percentage Mode
NEXT	Press to move between menu selections

Table 3-9. Statistics Mode F-Key Commands

3.11.2 Selecting Statistic Mode

1. Press , then press **F4**. The current operation mode displays.
2. If needed, press **F4**. The currently displayed mode flashes.
3. Use **F1** or **F2** to scroll to **STAT**.
4. Press **F4**, then press  to return to operation mode.
5. Choose to clear all of the data by pressing **F3** (**YES**) or to retain the data by pressing **F4** (**NO**).
6. Place product on the weighing pan and press  to store the weight.
7. Remove the specimen.
8. To collect and store more weighing data, repeat Steps 6 and 7 until the required number of data items are collected.
9. To display the statistical operation results, press **F2** (**DATA**). The statistical operation data menu displays.
10. Press **F1** or **F2** to toggle through the menu items available. Press **F4** (**RET**) when ready to continue.

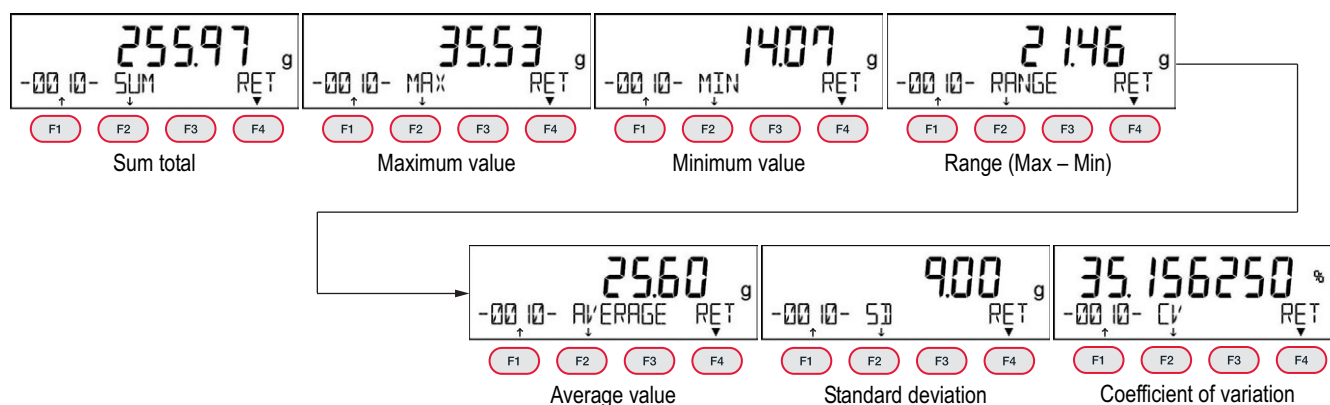


Figure 3-9. Statistics Mode – Data Menus

3.12 Animal Mode

The balance can accurately weigh animals and other live samples that can move by setting weight variations. When the live sample and movement fit within the set value range, the weight is held and displays. This mode is not available for Legal for Trade balances.

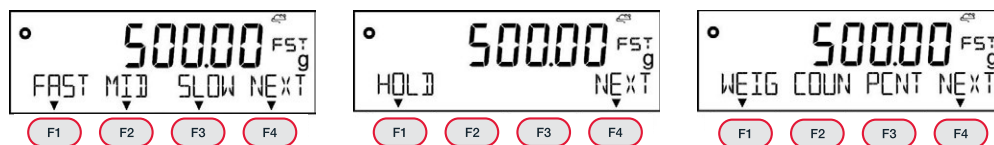





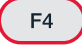







Figure 3-10. Animal Mode F-Keys

Command	Description
FAST	Used for animals that continue to move quickly while on the balance
MID	Used for animals that continue some movements while on the balance
SLOW	Used for animals that are calm while on the balance
HOLD	Manually press to hold the weight
WEIG	Press to display Weigh Mode
COUN	Press to display Count Mode
PCNT	Press to display Percentage Mode
NEXT	Press to move between menu selections

Table 3-10. Animal Mode F-Key Commands

When the external output is activated, the output condition is fixed as following:



- Output once after the indication is held except when the **HOLD** is pushed
- Output once after  is pressed during the indication is held

1. Press , then press . The current operation mode displays.
2. If needed, press . The currently displayed mode flashes.
3. Use  or  to scroll to **STAT**.
4. Press , then press  to return to operation mode.
5. Select the activity level by pressing  (**FAST**),  (**MID**) or  (**SLOW**).
6. Place the animal on the weighing pan. After the weight stabilizes within the set range the weight and **HLD** displays.
7. Remove the animal and the tare is automatically subtracted.



NOTE: For NTEP/Measurement Canada models, the zero-point automatically adjusts after the animal is removed.

3.12.1 Weigh Animal Using Manual Hold

1. Press  to display the **HOLD** menu.
2. Place the animal on the balance.
3. Press  to hold the weight reading at the time **HOLD** was pressed. The weight and **HLD** displays.
4. Remove the animal and the tare is automatically subtracted.



NOTE: For NTEP/Measurement Canada models, the zero-point automatically adjusts after the animal is removed.

3.13 Formulation Mode



NOTE: The formulation mode is only available for NTEP only models.

NTEP/Measurement Canada models do not support formulation mode.

Formulation Mode stores and refers the weight of each component compounded. This mode is not available for Legal for Trade balances.

- Only mg or g can be used
- Up to 30 components can be stored
- Preset tare function cannot be used

The output timing is set to once at stable or once immediately after is pressed, regardless of the setting value of External input/output function in the Condition Menu.

Settings:

- ON – Once at stable after key is pushed
- OFF – Once immediately after key is pushed



Figure 3-11. Formulation Mode F-Keys

Command	Description
END	Press to end formulation
DATA	Display data menu options
WEIG	Press to display Weigh Mode
COUN	Press to display Count Mode
PCNT	Press to display Percentage Mode
NEXT	Press to move between menu selections

Table 3-11. Formulation Mode F-Key Commands

1. Press , then press **F4** . The current operation mode displays.
2. If needed, press **F4** . The currently displayed mode flashes.
3. Use **F1** or **F2** to scroll to **FORM**.
4. Press **F4** , then press to return to operation mode. **MEM CLEAR YES NO** displays.
5. Choose to clear the memory data by pressing **F3** (**YES**) or save it by pressing **F4** (**NO**).
6. Place the weighing container on the balance and press to store the tare weight.
7. Put a specimen on the container and press to store the sample. **Now Sampling** flashing, then **Tare** followed by Net IMPORTED and the number of specimens.
8. Repeat steps 6-7 for all samples to be compounded.

Clear Data from Operation Mode

To end the formulation process and remove the stored data:




1. Press **F1** to end the formulation. **MEM CLEAR YES NO** displays.
2. Press **F3** (**YES**) to clear all data.

3.13.1 Data Review

1. To review the data, press **F2**.
2. Press **F1** or **F2** to scroll the data for listed formula.
 - Net – Net weight for listed formula
 - Tare – Tare for listed formula
 - Net Total – Total net weight of all net weights
 - Tare Total – Total tare weight of all net weights
3. Press **F3** to change the formula.
4. Press **F4** to return to operation mode.

3.14 Unit Setting

There are many units types available. In Legal for Trade balances, only **g** and **ct** are available. See [Section 4.3.2 on page 44](#) for more information and a complete list of available units.

1. Press .  displays.
2. Press **F4**. **11 MODE** and the current operation mode displays.
3. Press **F1** or **F2** to scroll to **UNIT**.
4. Press **F4** to change the display unit.
5. Press **F1** or **F2** to scroll to the desired unit. See [Table 4-2 on page 44](#).
6. Press **F4** to save the operation mode.
7. Press  to return to operation display.

3.15 Comparator Function

The Comparator Function is used to preset threshold values (limits) and determine if a measured value is within the preset range. This function can be used in Weigh Mode, Percentage Mode, Counting Mode and Multiplied by Coefficient Mode.

How to Perform Discrimination

Set the lower and the upper limits and whether the weight of a specimen to be weighed is **LOW** (lower than the lower limit), **OK** (appropriate) or **HIGH** (higher than the upper limit).

For setup information, refer to [Section 4.3.4 on page 46](#).

Discrimination	Single point setting (lower limit)	Single point setting (upper limit)	Two-point setting (upper and lower limits)
Over the upper limit	< OK > Blinking	< HIGH > Blinking	< HIGH > Blinking
Appropriate amount	< OK > Blinking	< OK > Blinking	< OK > Blinking
Below the lower limit	< LOW > Blinking	< OK > Blinking	< LOW > Blinking

Table 3-12. Messages

The discrimination is performed according to the following criteria:

Absolute value: The discrimination is performed based on the upper and lower limit values that have been set in advance.

Relative value: A reference numeric value is set in advance, and the discrimination is performed based on the range defined by the upper and lower limit values that have been set for the reference numeric value.

Example:

- Two-point (upper and lower limits) setting
- Reference value = 1000.00g
- Lower limit value = 900.00 g, Upper limit value = 1200.00 g

Discrimination Method	Reference Value	Lower Limit Value	Upper Limit Value
	1000.00 g	900.00 g	1200.00 g
Absolute value		900.00 g	1200.00 g
Relative value	1000.00 g	-100.00 g	200.00 g










Table 3-13. Example


3.16 Addition Function

The Addition Function is used to weigh several specimens in sequence and indicate the total value. The addition function can be used in Weighing Mode, Percentage Mode, Counting Mode and Multiplied by Coefficient Mode.

The addition function includes two ways of calculating method.

- Addition accumulating: specimens are weighed, removed and new specimens weighed
- Net adding function: specimens are weighed, more specimens added without removing previous

1. Press .  displays.
2. Press . **11 MODE** and the current operation mode displays.
3. Press  or  to scroll to **ADDITION**.
4. Press  to enter the menu. **ACTIVATE** displays with current setting.
5. Press . The current setting begins to flash.
6. Press  or  to select desired setting.

7. Press **F4** to save the setting.
8. Press **F1** or **F2** to scroll to **OPERATION**.
9. Press **F4**. The current setting begins to flash.
10. Press **F1** or **F2** to select desired setting.
11. Press **F4** to save the setting.
12. Press **F1** or **F2** to scroll to **DIRECTION**.
13. Press **F4**. The current setting begins to flash.
14. Press **F1** or **F2** to select desired setting.
15. Press **F4** to save the setting.
16. Press  to return to operation display.

3.16.1 Weighing with Plus Side Addition

Set ADD to **F1** and TOTL to **F2** for this function. See [Section 4.8.2 on page 59](#).

1. Place a specimen to be weighed on the balance. Wait for ***** to display.
2. Press **F1**. The weighed value is stored and **SUM TOTAL** displays momentarily.

For Addition Accumulating

1. Remove the first specimen and wait for **0.00** to display.
2. Place the next specimen on the balance. Wait for ***** to display.
3. Press **F1**. The weighed value is stored and **SUM TOTAL** displays momentarily.
4. Repeat this procedure until all additions have been completed.

For Net Addition

1. Perform steps 1-2 under [Section 3.16.1](#).
2. Without removing the previous specimen, add the next specimen to be weighed.
3. Once ***** displays, press **F1**. **SUM TOTAL** displays.
4. The weight is displayed, followed by an automatic tare. Repeat this procedure until all additions have been completed.




NOTE: For NTEP/Measurement Canada models, the weight is displayed, followed by the automatic zero-point adjustment.

5. Press **F2** (**TOTL**) to indicate the total value.
6. Press **F3** (**DEL**) to delete the total value.

3.16.2 Weighing with Minus Side Addition

Set ADD to **F1** and TOTL to **F2** for this function. See [Section 4.8.2 on page 59](#).

1. Place a first specimen to be weighed, ***** displays.
2. Press  **T**, **0.00** displays.



NOTE: Press  **0** for NTEP/Measurement Canada models.

For Addition Accumulating

1. Perform steps 1-2 under [Section 3.16.2](#).
2. Remove the specimen, ***** displays.
3. Press **F1**. The weighed value is stored and **SUM TOTAL** displays momentarily.



NOTE: Press  **0** for NTEP/Measurement Canada models after completing step 3.

4. Repeat this procedure until all additions have been completed.

For Net Addition

1. Perform steps 1-2 under [Section 3.16.2](#).
2. Add a specimen to be weighed.
3. Once ***** displays, press **F1**. **SUM TOTAL** displays momentarily.
4. The weight is displayed, followed by an automatic tare. Repeat this procedure until all additions have been completed.



NOTE: For NTEP/Measurement Canada models, the weight is displayed, followed by the automatic zero-point adjustment.

3.16.3 View/Delete Total Value

Press **F2** (**TOTL**) to indicate the total value.

Press **F3** (**DEL**) to delete the total value.

3.17 Tare-Subtraction Reminder Function



NOTE: The Tare-Subtraction Reminder Function is only available for NTEP only models.



If the Tare-Subtraction Reminder is activated, **PUSH TARE** displays when a container is placed on the balance.

If the Zero-Point-Adjustment Reminder is activated at the same time, it has priority.

There are two modes in the Tare-Subtraction Reminder function:

- **T Reminder 1** indicates the weight display is over the zero-point-adjustment range
- **T Reminder 2** indicates the weight display is over the zero-point-adjustment range before tare subtraction, and the net display is negative after tare subtraction

3.17.1 Set Tare-Subtraction Reminder Function

1. Press  **Menu**,  displays.
2. Press **F4**. **11 MODE** and the current operation mode displays.
3. Press **F1** or **F2** to scroll to **T REMINDER**.

4. Press **F4**. The current setting begins to flash.
5. Press **F1** or **F2** to select desired setting.
6. Press **F4** to save the setting.
7. Press **Menu** to return to operation display.

3.18 Zero-Point-Adjustment Reminder Function




NOTE: The Zero-Point-Adjustment Reminder Function is only available for NTEP only models.

If the Zero-Point-Adjustment Reminder is activated, **PUSH ZERO** is displayed if the load returns to within the zero point adjustment range after the load was over the range.

Place the specimens on the weighing pan. When they are removed, PUSH ZERO displays. Press **Zero** and **0** displays

3.18.1 Set Zero-Point-Adjustment Reminder Function


1. Press **Menu**.  displays.
2. Press **F4**. **11 MODE** and the current operation mode displays.
3. Press **F1** or **F2** to scroll to **Z REMINDER**.
4. Press **F4**. The current setting begins to flash.
5. Press **F1** or **F2** to select desired setting.
6. Press **F4** to save the setting.
7. Press **Menu** to return to operation display.

3.19 Stabilization Wait Function

The Stabilization Wait Function indicates when the weighed value displays after a zero-point adjustment or tare (only for NTEP only models), either after or before the weighed value stabilizes.











- OFF: function is not available
- ON: balance always waits for stabilization before displaying weighed value after the zero-point adjustment or tare

3.19.1 Set Stabilization Wait Function

1. Press **Menu**.  displays.
2. Press **F4**. **11 MODE** and the current operation mode displays.
3. Press **F1** or **F2** to scroll to **WT STABLE**.
4. Press **F4**. The current setting begins to flash.
5. Press **F1** or **F2** to select desired setting.
6. Press **F4** to save the setting.
7. Press **Menu** to return to operation display.

3.20 Bar Graph Display

The bar graph is displayed above the weight display when set to on.











1. Press .  displays.
2. Press . **11 MODE** and the current operation mode displays.
3. Press  or  to scroll to **BARGRAPH**.
4. Press . The current setting begins to flash.
5. Press  or  to select **ON** or **OFF**.
6. Press  to save the setting.
7. Press  to return to operation display.

3.21 Back Light Display

The back light can be set to **3MIN**, **5MIN**, **10MIN**, **30MIN**, **ON** or **OFF**. The back light will be on for the set number of minutes or always on if **ON** is selected. If **OFF** is selected, there is no back light.

This function does not work under if:

- a menu is displayed
- a specimen is placed on the weighing pan and the display is not stable











1. Press .  displays.
2. Press . **11 MODE** and the current operation mode displays.
3. Press  or  to scroll to **BACKLIGHT**.
4. Press . The current setting begins to flash.
5. Press  or  to select desired time.
6. Press  to save the setting.
7. Press  to return to operation display.

3.22 Auto Off Function

The Auto Off function is used to set an amount of time for the balance to remain on with no activity. Once the set time is reached the balance will automatically turn off. It can be set to **3MIN**, **5MIN**, **10MIN**, **30MIN** or **OFF**. If set to **OFF**, the balance remains on indefinitely.

This function does not work under if:











- a menu is displayed
- a specimen is on the weigh pan and the display is not stable

1. Press .  displays.
2. Press . **11 MODE** and the current operation mode displays.
3. Press  or  to scroll to **BACKLIGHT**.
4. Press . The current setting begins to flash.
5. Press  or  to select desired time.
6. Press  to save the setting.
7. Press  to return to operation display.

3.23 Simple Self Counting System (SCS) Method

1. Put a set number of samples in place.
2. Put up to two times the set number of additional samples in place. The balance will automatically update the average sample weight.
3. Repeating this allows accurate counting.





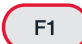







3.23.1 Set SCS Method

1. Press .  displays.
2. Press . **11 MODE** and the current operation mode displays.
3. Press  or  to scroll to **SIMPLE SCS**.
4. Press . The current setting begins to flash.
5. Press  or  to select ON or OFF
6. Press  to save the setting.
7. Press  to return to operation display.

4.0 Configuration

4.1 General Navigation

Use the keys to navigate through the menus and settings.

- Press  to enter menu structure.
- Press  or  to scroll through the main menus.
- Press  to enter a displayed menu.
- Press  or  to scroll through settings.
- Press  to enter the displayed setting. The current selection will flash.
- Press  or  to scroll through selections.
- Press  to select the displayed selection, it stops flashing.
- Press  to return to main menu.
- Press  at anytime to save changes and return to operation mode.




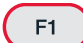

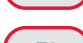

NOTE: The functions of the F-keys \uparrow , \downarrow , \leftarrow , \rightarrow , \leftarrow , or \blacktriangledown appear on the display directly above the corresponding F-key. Shortcuts for various modes/functions can be assigned to F-keys ([Section 4.8 on page 58](#)).

4.2 Numeric Value Entry





Figure 4-1. Front Panel Display


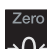
4.2.1 General Entry for All Models

-  Enter menu, cancel an input value, save parameter changes and return to previous menu, return to operation mode
-  = \uparrow Increments numeric values. Scrolls through menus/selections
-  = \downarrow Decrements numeric value. Scrolls through menus/selections
-  = \leftarrow Select digits and returns to previous level
-  = \rightarrow or \blacktriangledown Enters the value or a level

4.2.2 For NTEP Only Models

-  Use for changing polarity
-  Input a decimal point in **Multiplied by Coefficient** mode and **Specific Gravity** mode

4.2.3 For NTEP / Measurement Canada Models

-  Use for changing polarity
-  Input a decimal point in **Specific Gravity** mode
*Zero key on the right side of the display only

4.3 Applications Menu

The Applications Menu is used to set operational parameters.

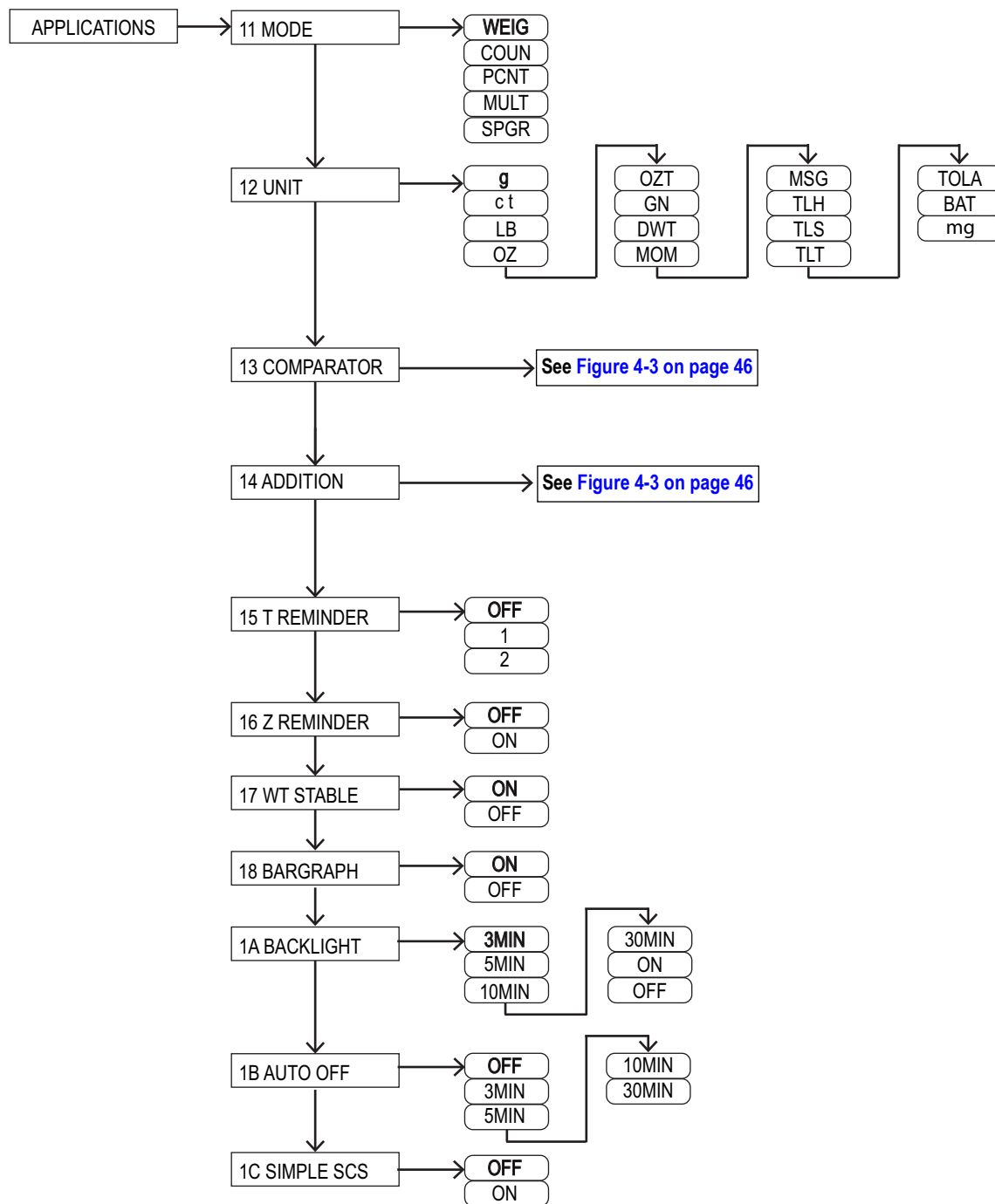


Figure 4-2. Applications Menu Layout

Menu	Parameters	Description
11 MODE	Select operation mode, see Section 4.3.1 on page 43	
	WEIG	Weigh Mode
	COUN	Counting Mode
	PCNT	Percentage Mode
	MULT*	Multiplied by a Coefficient
	SPGR	Specific Gravity Mode
	STAT*	Statistics Mode
	ANIM*	Animal Weighing Mode
	FORM	Formulation Mode; only for NTEP only models
12 UNIT***	Set the unit for weight, see Section 4.3.2 on page 44	
	g	Gram
	c t	Carat
	LB**	Pound
	OZ	Ounce
	OZT*	Troy Ounce
	GN**	Grain
	DWT*	Pennyweight
	MOM*	Momme
	MSG*	Mesghal
	TLH*	Tael (Hong Kong)
	TLS*	Tael (Singapore, Malaysia)
	TLT*	Tael (Taiwan)
	TOLA*	Tola
	BAT*	Baht
	mg*	Milligram
13 COMPARATOR	Comparator Function; see Figure 4-3 on page 46	
14 ADDITION*	Adding Function; see Section 4.3.5 on page 47	

Menu	Parameters	Description
15 T REMINDER	Tare on Reminder; only for NTEP only models; see Section 3.17 on page 37	
	OFF	Disabled
	1	Activate the MODE 1
	2	Activate the MODE 2
16 Z REMINDER	Zero on Reminder; only for NTEP only models; see Section 3.18 on page 38	
	OFF	Disabled
	ON	Enabled
17 WT STABLE*	Stability Waiting; see Section 3.19 on page 38	
	ON	Enabled
	OFF	Disabled
18 BARGRAPH	Bar Graph Indication; see Section 3.20 on page 39	
	ON	Enabled
	OFF	Disabled
1A BACKLIGHT	Back light on timer, see Section 3.21 on page 39	
	3MIN	3 minutes
	5MIN	5 minutes
	10MIN	10 minutes
	30MIN	30 minutes
	ON	Always on
	OFF	Disabled
1A AUTO OFF	Auto Power Off, time the balance is inactive before it powers off; see Section 3.22 on page 40	
	OFF	Disabled
	3MIN	3 minutes
	5MIN	5 minutes
	10MIN	10 minutes
	30MIN	30 minutes
	30MIN	30 minutes
1C SIMPLE SCS	Simplified SCS; see Section 3.23 on page 40	
	OFF	Disabled
	ON	Enabled

Table 4-1. Applications Menu Parameters











NOTE: *Parameter menu/options not available on Legal for Trade verified NTEP/Measurement Canada models.

**Pound (LB) and Grain (GN) are not selectable on Legal for Trade verified TE-1501NC models.

***12 UNIT is not selectable for Legal for Trade verified TE-8200NC models and only Gram (g) is available.

4.3.1 Operation Modes

Set the balance to desired operation mode.

1. Press .  displays.
2. Press . **11 MODE** and the current operation mode displays.
3. Press  to change the mode if needed.
4. Press  or  to scroll to the desired mode of operation.
5. Press  to save the operation mode.
6. Press  to return to operation display.

4.3.2 NTEP Only Model Units

Use Table 4-2 when setting the display units. Not all units are available in Legal for Trade balances.






NOTE: For NTEP only models, $e \neq d$. The balance's readability (d) is 10x the NTEP verified readability (e). (e.g. The TE-223 has a readability (d) of 0.001 g and an NTEP verified readability (e) of 0.01 g.)

Display	Unit	Conversion Coefficient	Model Capacity and Readability by Unit				
			TE-223	TE-623	TE-3202	TE-6202	TE-15001
g	gram	1.00000000E+00	220 0.001	620 0.001	3200 0.01	6200 0.01	15000 0.1
ct	carat	5.00000000E+00	1100 0.01	3100 0.01	16000 0.1	31000 0.1	75000 1
lb	pound	2.20462260E-03	0.48 0.00001	1.3 0.00001	7 0.0001	13 0.0001	33 0.001
oz	ounce	3.52739610E-02	7.7 0.0001	21 0.0001	110 0.001	210 0.001	520 0.01
ozt	troy ounce	3.21507460E-02	7 0.0001	19 0.0001	100 0.001	190 0.001	480 0.01
GN	grain	1.54323580E+01	3300 0.1	9500 0.1	49000 1	95000 1	230000 10
dwt	pennyweight	6.43014930E-01	140 0.001	390 0.001	2000 0.01	3900 0.01	9600 0.1
mom	momme	2.66666670E-01	58 0.001	160 0.001	850 0.01	1600 0.01	4000 0.1
MSG	mesghal	2.16999761E-01	47 0.001	130 0.001	690 0.01	1300 0.01	3200 0.1
t:H	Hong Kong tael	2.67172510E-02	5.8 0.0001	16 0.0001	85 0.001	160 0.001	400 0.01
t:S	Singapore, Malaysia tael	2.64554710E-02	5.8 0.0001	16 0.0001	84 0.001	160 0.001	390 0.01
t:T	Taiwan tael	2.66666670E-02	5.8 0.0001	16 0.0001	85 0.001	160 0.001	400 0.01
to	tola	8.57353240E-02	18 0.0001	53 0.0001	270 0.001	530 0.001	1200 0.01
BA	baht	6.59630607E-02	14 0.0001	40 0.0001	210 0.001	400 0.001	980 0.01
mg	milligram	1.00000000E+03	220000 1	620000 1	3200000 10	6200000 10	15000000 100

Table 4-2. Unit Conversion Weighing Capacity and Readability (NTEP Only Models)

To set the desired unit displayed:

1. Press .  displays.
2. Press **F4**. **11 MODE** and the current operation mode displays.
3. Press **F1** or **F2** to scroll to **UNIT**.
4. Press **F4** to change the display unit.
5. Press **F1** or **F2** to scroll to the desired unit.
6. Press **F4** to save the operation mode.
7. Press  to return to operation display.

4.3.3 NTEP/Measurement Canada Model Units

Use [Table 4-3](#) when setting the display units. Not all units are available in Legal for Trade balances. See procedure steps on the previous page to set the desired unit displayed.



NOTE: For NTEP/Measurement Canada models, $e = d$. The balance's readability (d) is equal to the verified readability (e). (e.g. The TE-322NC has a readability (d) of 0.01 g and a verified readability (e) of 0.01 g.)

Display	Unit	Conversion Coefficient	Model Capacity and Readability by Unit		
			TE-322NC	TE-1501NC	TE-8200NC
g	gram***	1.00000000E+00	320 0.01	1500 0.1	8200 1
ct	carat	5.00000000E+00	1600 0.1	7500 1	41000 10
lb	pound**	2.20462260E-03	0.7 0.0001	3.3 0.001	18 0.01
oz	ounce	3.52739610E-02	11 0.001	52 0.01	280 0.1
oz t	troy ounce*	3.21507460E-02	10 0.001	48 0.01	260 0.1
GN	grain**	1.54323580E+01	4900 1	23000 10	120000 100
dwt	pennyweight*	6.43014930E-01	200 0.01	960 0.1	5200 1
mom	momme*	2.66666670E-01	85 0.01	400 0.1	2100 1
MSG	mesghal*	2.16999761E-01	69 0.01	320 0.1	1700 1
t : H	Hong Kong tael*	2.67172510E-02	8.5 0.001	40 0.01	210 0.1
t : S	Singapore, Malaysia tael*	2.64554710E-02	8.4 0.001	39 0.01	210 0.1
t : T	Taiwan tael*	2.66666670E-02	8.5 0.001	40 0.01	210 0.1
to	tola*	8.57353240E-02	27 0.001	120 0.01	700 0.1
BA t	baht*	6.59630607E-02	21 0.001	98 0.01	540 0.1
mg	milligram*	1.00000000E+03	320000 10	1500000 100	8200000 1000

Table 4-3. Unit Conversion Weighing Capacity and Readability (NTEP/Measurement Canada Models)



NOTE: *Unit not available on Legal for Trade verified NTEP/Measurement Canada models.

**Pound (LB) and Grain (GN) are not selectable on Legal for Trade verified TE-1501NC models.

***Gram (g) is the only available unit for Legal for Trade verified TE-8200NC models.

4.3.4 Comparator Menu

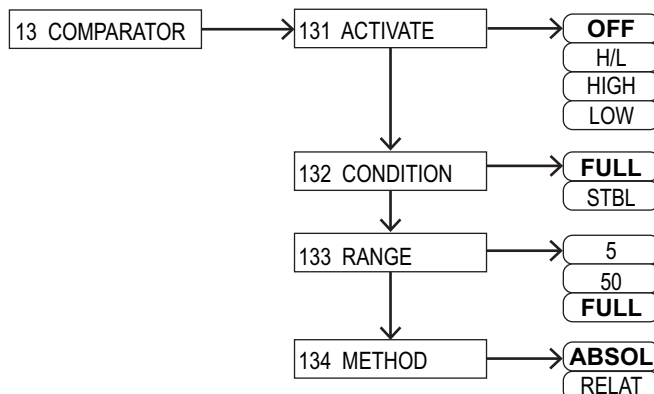
















Figure 4-3. Applications Comparator Menu Layout

Menu	Parameters	Description
131 ACTIVATE	Activates the comparator function	
	OFF	Disabled
	H/L	Valid upper and lower limits
	HIGH	Valid upper limit
	LOW	Valid lower limit
132 CONDITION	Discriminant condition	
	FULL	Always
	STBL	Only at stable times
133 RANGE	Discriminant range	
	5	+5 (e/d) or more
	50	+50 (e/d) or more
	FULL	Entire area
134 METHOD	Discriminant method	
	ABSOL	Absolute value method
	RELAT	Relative value method

Table 4-4. Applications Comparator Parameters

See [Section 3.15 on page 35](#) for user information.

1. Press .  displays.
2. Press . **11 MODE** and the current operation mode displays.
3. Press  or  to scroll to **COMPARATOR**.
4. Press  to enter the menu. **ACTIVATE** displays with current setting.
5. Press . The current setting begins to flash.
6. Press  or  to select desired setting (**OFF**, **H/L**, **HIGH** or **LOW**).
7. Press  to save the setting.
8. Press  or  to scroll to **CONDITION**.
9. Press . The current setting begins to flash.

10. Press **F1** or **F2** to select desired setting.
11. Press **F4** to save the setting.
12. Press **F1** or **F2** to scroll to **RANGE**.
13. Press **F4**. The current setting begins to flash.
14. Press **F1** or **F2** to select desired setting.
15. Press **F4** to save the setting.
16. Press **F1** or **F2** to scroll to **METHOD**.
17. Press **F4**. The current setting begins to flash.
18. Press **F1** or **F2** to select desired setting.
19. Press **F4** to save the setting.
20. Press  to return to operation display.

4.3.5 Addition

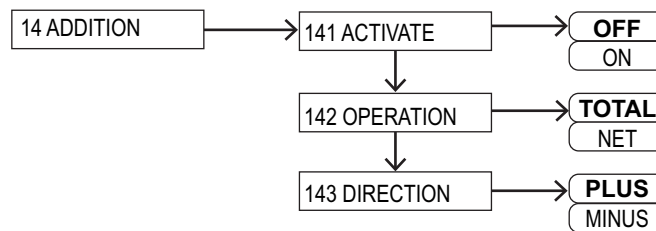




Figure 4-4. Applications Addition Menu Layout




NOTE: 14 ADDITION is not selectable for Legal for Trade verified NTEP/Measurement Canada models.

Menu	Parameters	Description
131 ACTIVATE	Activates the addition function	
	OFF	Disabled
	ON	Enabled
132 OPERATION	Adding operation	
	TOTAL	Addition accumulated
	NET	Net addition
133 DIRECTION	Adding direction	
	PLUS	Plus side addition
	MINUS	Minus side addition

Table 4-5. Applications Addition Parameters

1. Press .  displays.
2. Press **F4**. **11 MODE** and the current operation mode displays.
3. Press **F1** or **F2** to scroll to **ADDITION**.

4. Press **F4** to enter the menu. **ACTIVATE** displays with current setting.
5. Press **F4**. The current setting begins to flash.
6. Press **F1** or **F2** to select desired setting.
7. Press **F4** to save the setting.
8. Press **F1** or **F2** to scroll to **OPERATION**.
9. Press **F4**. The current setting begins to flash.
10. Press **F1** or **F2** to select desired setting.
11. Press **F4** to save the setting.
12. Press **F1** or **F2** to scroll to **DIRECTION**.
13. Press **F4**. The current setting begins to flash.
14. Press **F1** or **F2** to select desired setting.
15. Press **F4** to save the setting.
16. Press  to return to operation display.

4.4 Performance Menu

Set the balance display stability, response and Zero Track speed.

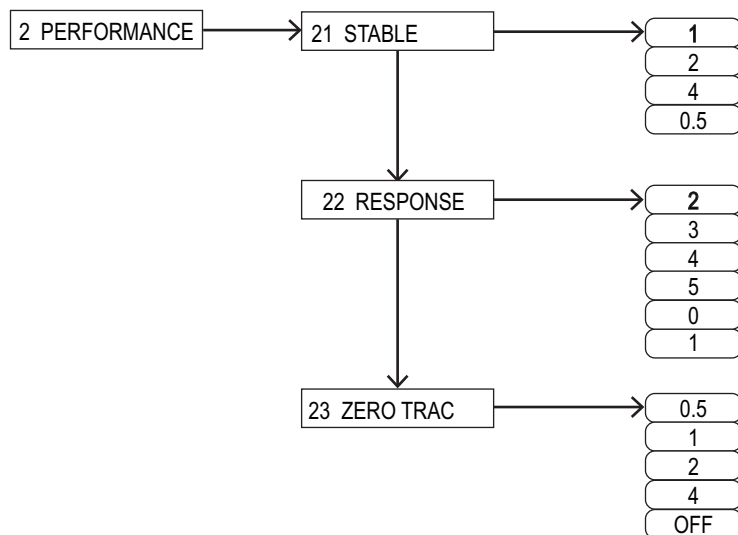


Figure 4-5. Performance Menu Layout

Menu	Parameters	Description
21 STABLE	Stability discrimination width	
	1*	1d
	2	2d
	4	4d
	0.5	0.5d
22 RESPONSE	Response speed	
	0*	Sensitive Mode
	1	Fast
	2	Medium fast
	3	Medium
	4	Medium slow
	5	Slow
23 ZERO TRAC	Zero Tracking	
	0.5*	0.5d
	1	1d
	2	2d
	4	4d
	OFF	Invalid

*Fixed setting for NTEP/Measurement Canada models.

Table 4-6. Performance Menu Parameters











4.4.1 Stability

When the larger numeric value is set in this setting menu, the laser stability judgment is applied and the balance displays **O** in more unstable conditions.



NOTE: Selections 2 and 4 cannot be selected in Legal for Trade balances.

For verified NTEP / Measurement Canada models, this function is not available and 21 STABLE is fixed at 1.

1. Press .  displays.
2. Press  or  to scroll to **PERFORMANCE**.
3. Press . **21 STABLE** and the current setting displays.
4. Press . The current setting begins to flash.
5. Press  or  to select desired setting.
6. Press  to save the setting.
7. Press  to return to operation display.









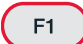



4.4.2 Response Speed

The larger the value is set in this menu, the more stable the balance display becomes in unstable conditions.



NOTE: Selections 1, 2 and 4 cannot be selected in Legal for Trade balances.

For verified NTEP / Measurement Canada models, this function is not available and 22 RESPONSE is fixed at 0.

1. Press .  displays.
2. Press  or  to scroll to **PERFORMANCE**.
3. Press . **21 STABLE** and the current setting displays.
4. Press  or  to scroll to **22 RESPONSE**.
5. Press . The current setting begins to flash.
6. Press  or  to select desired setting.
7. Press  to save the setting.
8. Press  to return to operation display.









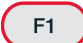



4.4.3 Zero Tracking

The Zero Tracking function makes it possible to automatically correct the zero-point fluctuation when **0** is displayed, through which the **0** display is maintained.



NOTE: Selections 1, 2 and 4 cannot be selected in Legal for Trade balances.

For verified NTEP / Measurement Canada models, this function is not available and 23 ZERO TRAC is fixed at 0.5.

1. Press .  displays.
2. Press  or  to scroll to **PERFORMANCE**.
3. Press . **21 STABLE** and the current setting displays.
4. Press  or  to scroll to **23 ZERO TRAC**.
5. Press . The current setting begins to flash.
6. Press  or  to select desired setting.
7. Press  to save the setting.
8. Press  to return to operation display.

4.5 User Information

Describes setting items related to the comparator function and preset tare weight.



NOTE: 31 PT MODE and 32 PT INPUT are only available in NTEP only models.
NTEP/Measurement Canada models do not support the tare function.

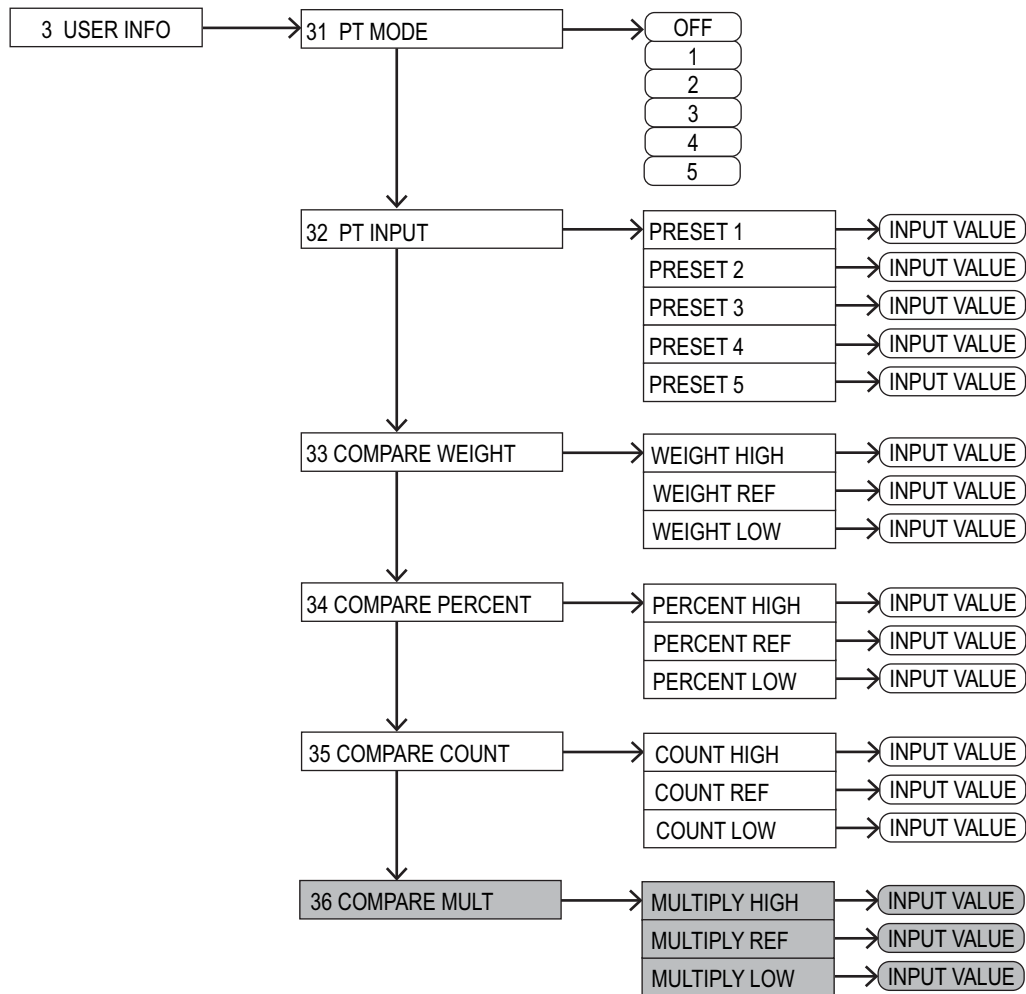


Figure 4-6. User Information Menu

Menu	Parameter	Description
31 PT MODE	Preset Tare Execution; only available for NTEP only models	
	Off	Invalid
	Setting 1-5	Execute
32 PT INPUT	Preset Tare Weight Setting; only available for NTEP only models	
	Preset 1-5	Setting value input
33 COMPARE WEIGHT	Weight Comparator	
	Weight High	Setting value input
	Weight Ref	
	Weight Low	

Menu	Parameter	Description
34 COMPARE PERCENT	% Comparator	
	Percent High	Setting value input
	Percent Ref	
	Percent Low	
35 COMPARE COUNT	Comparator Counting	
	Count High	Setting value input
	Count Ref	
	Count Low	
36 COMPARE MULT (Not available for verified Legal for Trade models)	Multiple Comparator	
	Multiply High	Setting value input
	Multiply Ref	
	Multiply Low	

Table 4-7. User Information Menu

4.5.1 Preset Tare Mode

A preset tare can be selected prior to weighing. See [Section 3.5.1 on page 24](#) for more information.















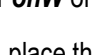





NOTE: The preset tare function is only available in NTEP only models.

4.5.2 Input Preset Tare Value

There are two ways of inputting a preset tare weight value

- Actual value setting method (onW): Weighing a sample with a scale and then making it a setting value
- Numeric value setting method (NUM): Inputting a setting value directly via key operation

To enter a preset tare:

1. Press .  displays.
2. Press  or  to scroll to **USER INFO**.
3. Press . **31 PT MODE** displays.
4. Press  or  to scroll to **32 PT INPUT**.
5. Press . **PRESET 1** displays.
6. Press  or  to select desired preset tare (1-5).
7. Press . **SET PRESET onW NUM** displays.
8. Select  for **onW** or  for **NUM**.
 If **onW** is selected, place the container on the scale, when weight is stable, press  to save
 If **NUM** is selected use  or  to enter known value, press  to save
9. Press  to return to operation display. **NET PT** display in the upper right.

4.5.3 Set the Discrimination Value of Comparator Function

There are two ways of inputting a reference value and upper and lower limit values

- Actual value setting method (onW): Weighing a sample with a scale and then making it a setting value
- Numeric value setting method (NUM): Inputting a setting value directly via key operation

The discrimination is performed according to the following criteria:

Absolute Value

The discrimination is performed based on the upper and lower limit values that have been set in advance.

Relative Value

A reference numeric value is set in advance, and the discrimination is performed based on the range defined by the upper and lower limit values that have been set for the reference numeric value.

Example:














Two-point (upper and lower limits) setting, Reference value = 1000.00 g

Lower limit value = 900.00 g, Upper limit value = 1200.00 g

Discrimination Method	Reference value	Lower limit value	Upper limit value
	1000.00 g	900.00 g	1200.00 g
Absolute value		900.00 g	1200.00 g
Relative value	1000.00 g	-100.00 g	200.00 g

Table 4-8. Relative Value Example

To enter a value in the comparator function:

1. Press .  displays.
2. Press  or  to scroll to **USER INFO**.
3. Press . **31 PT MODE** displays.
4. Press  or  to scroll to the comparator function to be set.
5. Press  to enter the function.
6. Press  or  to select the parameter to set.
7. Press . Current function and **onW NUM** display.
8. Select  for **onW** or  for **NUM**.

If **onW** is selected, place the container on the scale, when weight is stable, press  to save

If **NUM** is selected use  or  to enter known value, press  to save

9. Press  to return to operation display. **NET PT** display in the upper right.

Repeat this section for each of the comparator settings:

- Comparator setting for Weighing mode: **33 COMPARE WEIGHT**
- Comparator setting for Percentage mode: **34 COMPARE PERCENT**
- Comparator setting for Counting mode: **35 COMPARE COUNT**
- Comparator setting for Multiplied by Coefficient mode: **36 COMPARE MULT**

4.6 External Input/Output Functions

This function is used for communication through the external peripheral devices. There are RS-232C (D-SUB 9P) and USB (Type B) interface as standard equipment, and each interface slot for option. See [Section 6.0 on page 67](#) for setup.

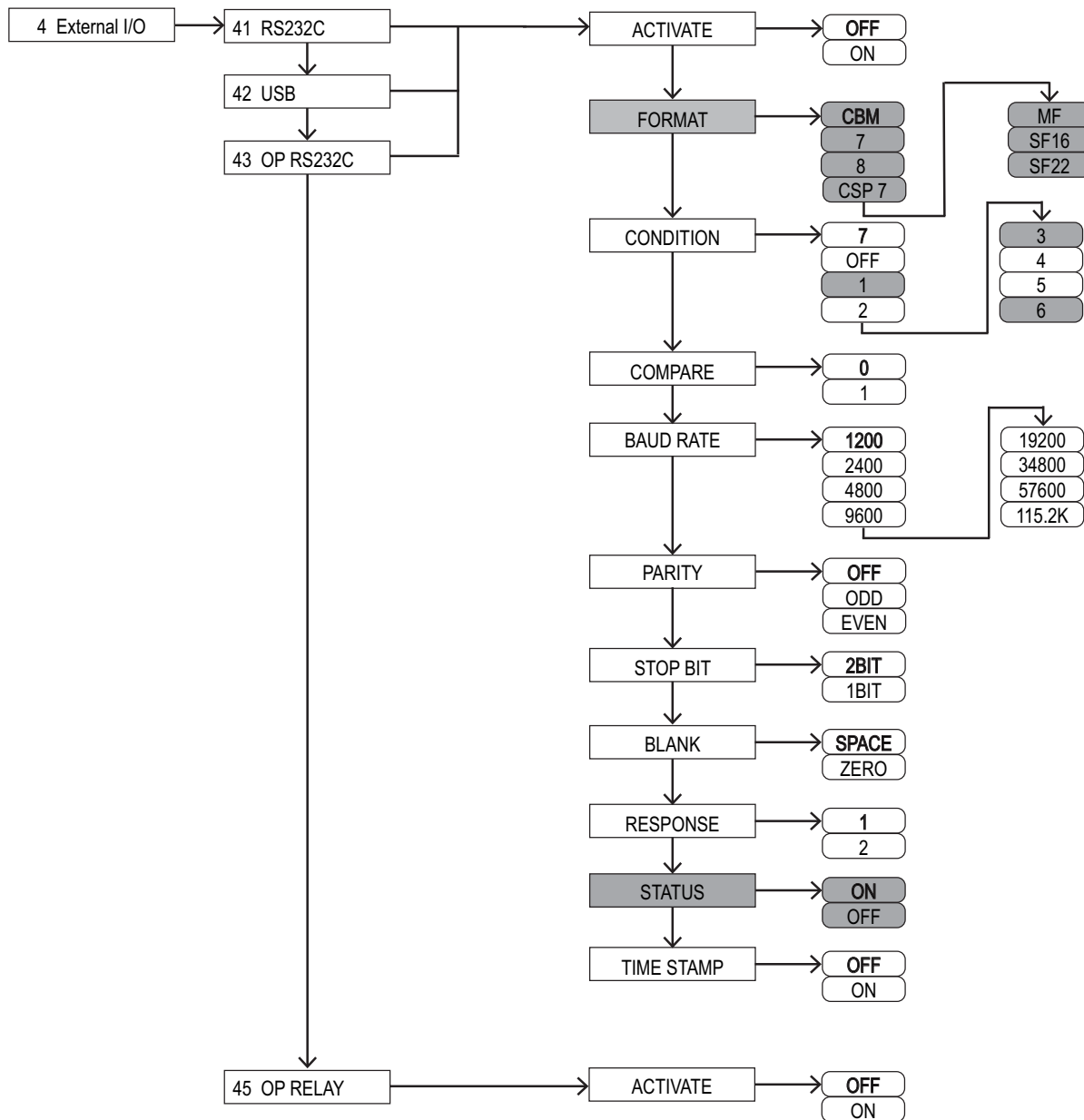


Figure 4-7. External Input/Output Menu Structure

Type	Menu	Parameters	Description
RS232C (Standard) USB (Standard) OP RS232C (Expanded Option)	ACTIVATE	OFF	Stop
		ON	Operation
	FORMAT	CBM	CBM Format
		6	6 Digit Format
		7	7 Digit Format
		8	8 Digit Format
		CSP 6	CSP Format 6 digits
		CSP 7	CSP Format 7 digits
	CONDITION	7	Push down (output key) for one time output at stable times
		OFF	Output Stop
		1	Continuous output at all times
		2	Continuous output as stable
		3	Push down (output key) for one time instant output
		4	Auto output
		5	One time output at stable times
		6	One time output at stable times
	COMPARE	0	Per output setting
		1	Output when discrimination result is okay or absent
	BAUD RATE	1200	1200 bps
		2400	2400 bps
		4800	4800 bps
		9600	9600 bps
		19200	19200 bps
		34800	34800 bps
		57600	57600 bps
		115.2K	115200 bps
	PARITY	OFF	NONE
		ODD	Odd Number
		EVEN	Even Number
	STOP BIT	2BIT	2 Bit
		1BIT	1 Bit
	BLANK	SPACE	Fill with a blank space (0 x 20)
		ZERO	Fill with 0 (0 x 30)
	RESPONSE	1	A00 Exx format
		2	ACK, NAK format
	STATUS	OFF	Not added
		ON	Append
OP LIMIT (Option) <i>Relay Output</i>	ACTIVATE	ON	Operation
		OFF	Stop

Table 4-9. External Input/Output Parameters

4.7 Lock Functions

Limitations can be imposed on key operation and in accessing menu items.

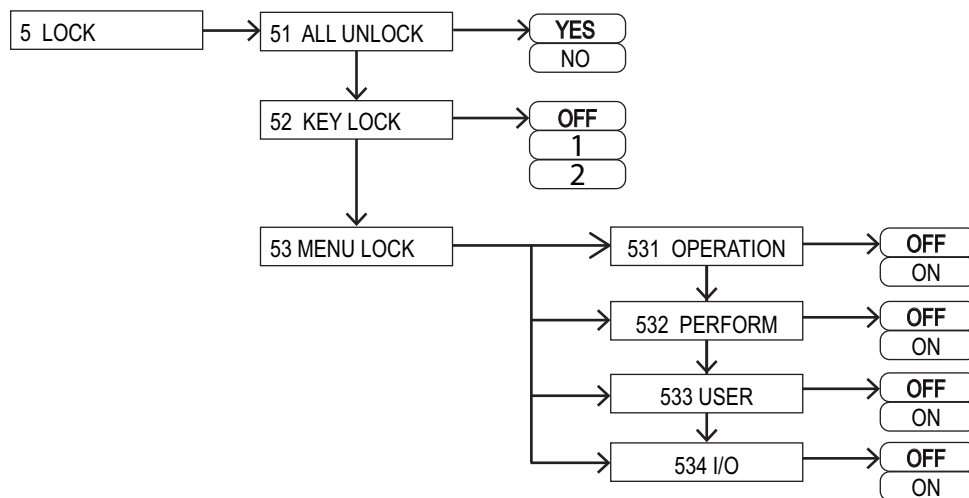


Figure 4-8. Lock Menu Structure


















4.7.1 Total Lock Release

Use the following steps to release all the locks that have been set.

1. Press . displays.
2. Press or until displays.
3. Press . displays.
4. Press . displays.
5. Press (YES) or (NO).
6. Press . To return to the operation mode.




















4.7.2 Key Lock Function

Use the following steps to lock key operation.

1. Press .  displays.
2. Press  or  until  displays.
3. Press .  displays.
4. Press  or  until  displays.
5. Press . **OFF** starts flashing.
6. Press  or  to desired setting.
 - **OFF**: all keys are available
 - **1**:  locked for power off
 - **2**: All keys are locked (except ; when in menu mode all keys are available)
7. Press . The chosen setting displays.
8. Press  to return to the operation mode.

4.7.3 Menu Lock Function

Various setting menus can be locked.

1. Press .  displays.
2. Press  or  until  displays.
3. Press .  displays.
4. Press  or  until  displays.
5. Press .  displays.
6. Press  or  to display the menus available for locking.
 - 531 OPERATION: Function related to the operation <1 APPLICATIONS>
 - 532 PERFORM: Function related to the performance <2 PERFORMANCE>
 - 533 USER: User information setting <3 USER INFO>
 - 534 I/O: External input/output functions <4 EXTERNAL I/O>
7. Press  at each menu to be locked/unlocked. The current setting flashes.
8. Press  or  to display **ON** or **OFF**.
9. Press . Setting stops flashing.
10. When all menus are set, press  to return to the operation mode.

4.8 Admin/Adjust Menu

Perform setting of the balance ID, the span adjustment and the date and time.

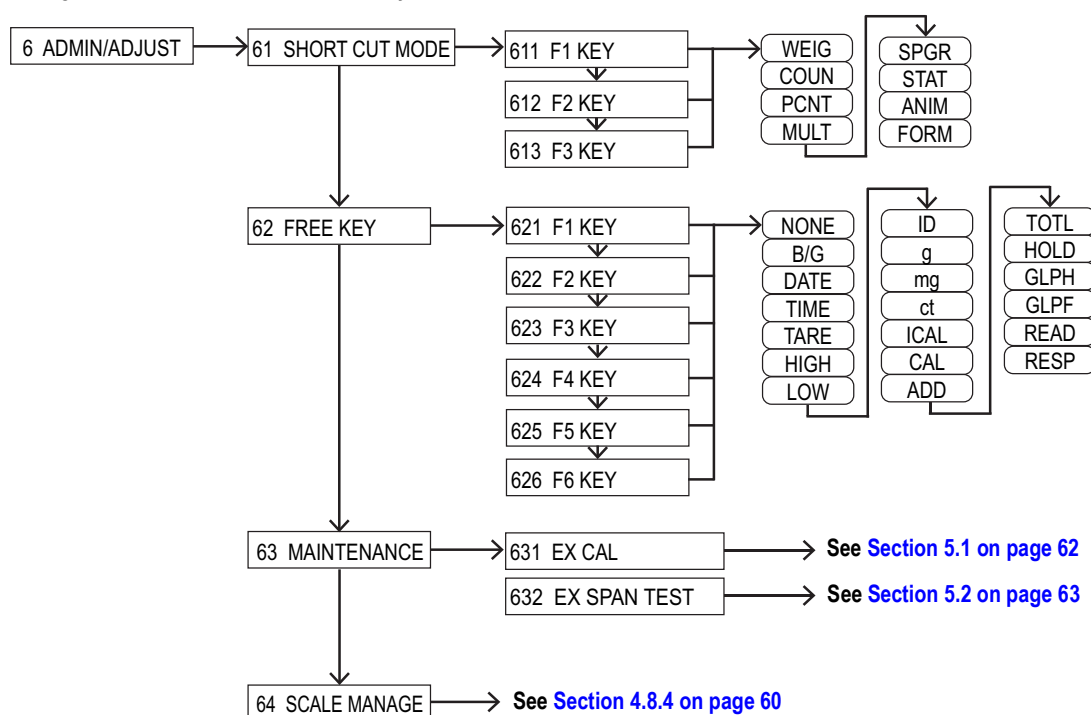


Figure 4-9. Controlling and Adjusting Menu Layout







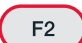






Menu	Sub-Menu	Parameters – Description
61 SHORT CUT MODE	Key assignment for mode selection	
	611 F1 KEY – WEIG (default)	WEIG – Weigh Mode; COUN – Counting Mode; PCNT – Percentage Mode
	612 F2 KEY – COUN (default)	MULT* – Multiplied by a Coefficient; SPGR – Specific Gravity Mode; STAT* – Statistics Mode
	613 F3 KEY – PCNT (default)	ANIM* – Animal Weighing Mode; FORM* – Formulation Mode (only in NTEP only models)
62 FREE KEY (NTEP only models)	Free key assignment; F1-F3 at layer 1; F4-F6 at layer 2;	
	621 F1 KEY – B/G (default)	B/G – brutto/gross; DATE – date indication; TIME – time indication; TARE – tare indication;
	622 F2 KEY – DATE (default)	HIGH – upper limit value; LOW – lower limit value; ID – ID number indication
	623 F3 KEY – TIME (default)	g – unit setting; mg* – unit setting; ct – unit setting; ICAL – internal span adjustment
	624 F4 KEY – TARE (default)	CAL* – external span adjustment; ADD – adding execute; TOTL – total indication
	625 F5 KEY – HIGH (default)	HOLD* – measurement indication hold; GLPH – GLP header printing
	625 F6 KEY – LOW (default)	GLPF – GLP footer printing; READ – designation of readability (d)
62 FREE KEY (NTEP/Measurement Canada models)	Free key assignment; F1-F3 at layer 1; F4-F6 at layer 2;	
	621 F1 KEY – DATE (default)	DATE – date indication; TIME – time indication; HIGH – upper limit value;
	622 F2 KEY – TIME (default)	LOW – lower limit value; ID – ID number indication; g – unit setting; mg* – unit setting;
	623 F3 KEY – ID (default)	ct – unit setting; CAL* – external span adjustment; ADD* – adding execute;
	624 F4 KEY – HIGH (default)	TOTL* – total indication; HOLD* – measurement indication hold;
	625 F5 KEY – LOW (default)	GLPH – GLP header printing; GLPF – GLP footer printing;
	625 F6 KEY – g (default)	RESP* – response speed; NONE – disabled
63 MAINTENANCE	Maintenance settings	
	631 EX CAL	External Calibration – Execute; see Section 5.1 on page 62
	632 EX SPAN TEST	External Span Calibration Test – Execute (rewrite); see Section 5.2 on page 63

Table 4-10. Controlling and Adjusting Parameters















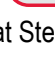


NOTE: *Parameters not available on Legal for Trade balances.

4.8.1 Set Short Cut Mode

1. Press .  displays.
2. Press  or  to scroll to **6 ADMIN/ADJUST**.
3. Press . The current sub-menu displays.
4. Press  or  to scroll to **61 SHORT CUT MODE**.
5. Press . The current F-Key displays.
6. Press . The current setting begins to flash.
7. Press  or  to select desired setting.
8. Press  to save the setting.
9. Press  to return to operation display.

4.8.2 Set Free Keys

1. Press .  displays.
2. Press  or  to scroll to **6 ADMIN/ADJUST**.
3. Press . The current sub-menu displays.
4. Press  or  to scroll to **62 FREE KEY**.
5. Press . The current F-Key displays.
6. Press  or  to scroll to desired F-Key.
7. Press . The current setting begins to flash.
8. Press  or  to select desired setting.
9. Press  to save the setting.
10. Repeat Steps 5-9 until all Free Keys are set.
11. Press  to return to operation display.

4.8.3 Maintenance Settings

The Maintenance Settings menu contains calibration parameters to perform an external calibration ([Section 5.1 on page 62](#)) or to perform an external span calibration test ([Section 5.2 on page 63](#)).

4.8.4 Balance Manage Menu

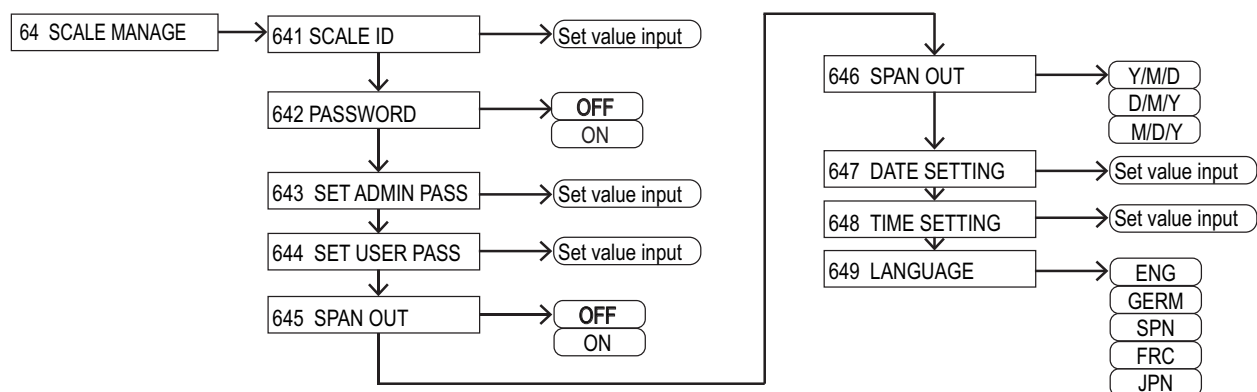


Figure 4-10. Controlling and Adjusting–Balance Manage Menu Layout

Menu	Parameters	Description
641 SCALE ID	Balance ID: Enter ID value	
642 PASSWORD	Password Control	
	OFF	Disabled
	ON	Enabled
643 SET ADMIN PASS	Administrators password registration–enter password value	
644 SET USER PASS	User password registration–enter password value	
645 SPAN OUT	Output of the span adjustment test	
	OFF	Disabled
	ON	Enabled
646 DISP DATE	Date display format	
	Y/M/D	Year, month, day
	D/M/Y	Day/month/year
	M/D/Y	Month/day/year
647 DATE SETTING	Enter date	
648 TIME SETTING	Enter time	
649 LANGUAGE	Printed Language	
	ENG	English
	GERM	Deutsch
	SPN	Espanol
	FRC	French
	JPN	Japanese
64A SPACING (Only available for NTEP only models)	Readability setting	
	1	1d
	2	2d
	5	5d
	10	10d
64B START CAL (Only available for NTEP only models)	Span adjustment with internal weight at power on	
	OFF	Disabled
	FORCE	Enabled
	SELEC	Selectable
64C DIRECT ST	Direct start setting	
	OFF	Disabled
	ON	Enabled
64D INITIALIZE	Initialize	
	YES	Cancel
	NO	Execute

Table 4-11. Controlling and Adjusting–Balance Manage Parameters

4.9 LCD Display Characters

4.9.1 7-Segment Characters

A	B	C	D	E	F	G	H	I	J	K	L	M	N	O
A	b	c	d	E	F	G	h	i	J	K	L	M	n	o
P	Q	R	S	T	U	V	W	X	Y	Z	c	comma	point	
P	q	r	s	t	u	v	w	x	y	z	c	,	.	
1	2	3	4	5	6	7	8	9	0	space	minus / hyphen			
1	2	3	4	5	6	7	8	9	0	-	-			

4.9.2 16-Segment Characters

A	B	C	D	E	F	G	H	I	J	K	L	M	N	O
A	B	C	D	E	F	G	H	I	J	K	L	M	N	O
P	Q	R	S	T	U	V	W	X	Y	Z				
P	Q	R	S	T	U	V	W	X	Y	Z				
b	c	d	g	l	m	n	o	t	w					
b	c	d	g	l	m	n	o	t	w					
1	2	3	4	5	6	7	8	9	0					
1	2	3	4	5	6	7	8	9	0					
asterisk	slash	left arrow	right arrow	space	plus	minus / hyphen								
*	/	←	→		+	-								
comma	point	percent	Degree Celsius											
,	.	%	°C											

5.0 Calibration

Calibration parameters are located in the Admin/Adjust Menu. Refer to [Section 4.8 on page 58](#) for menu location details.

A span adjustment ([Section 5.1](#)) is to “decrease” the difference between an indicated value and the true value (mass), and a span test ([Section 5.2 on page 63](#)) is to “check” the difference between an indicated value and the true value.

Regularly check the accuracy of the balance and perform calibrations as needed to maintain high-accuracy weighing. Since electronic balances are affected by the acceleration of gravity, an adjustment/test is needed at every weighing location. An adjustment/test is also needed when in use for a long period of time or if the accuracy indication no longer appears.

















IMPORTANT: Only use and calibrate the balance in areas free from environmental conditions that could affect the accuracy ([Section 1.3 on page 9](#)). Calibration weight must meet or exceed OIML F1 Tolerances.

5.1 Perform External Calibration

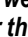

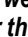
Doing a calibration and making a span adjustment can significantly affect weighing accuracy. Please read this procedure carefully before beginning.



NOTE: Calibration is only available before Legal for Trade verification.

1. Level balance by adjusting the feet, if needed ([Section 2.4 on page 15](#)).
2. Press .  displays.
3. Press  or  to scroll to **6 ADMIN/ADJUST**.
4. Press . **61 SHORT CUT MODE** displays.
5. Press  or  to scroll to **63 MAINTENANCE**.
6. Press . **631 EX CAL** displays.
7. Press . **CAL STARTED** displays, then **SELECT MIN** displays.
8. Press  or  to choose the division size (1, 2, 5 or 10).
9. Press . **SELECT WEIGHT** displays.
10. Press  or  to select the desired weight amount.





IMPORTANT: If a calibration weight less than the weighing capacity is used,   displays when weighing a sample that is more than two times heavier than the calibration weight.   indicates that the weighing accuracy is uncertain.



NOTE: There is the option to choose VAR to enter a custom weight amount.

Press  to select a digit,  or  to change a digit, and  to confirm amount.









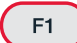


11. Press . **ON 0** displays, then **ON F** displays.
12. Place weight on pan. **ON F** flashes.
13. Once calibration is complete, **DATA SAVED** displays, then **631 EX CAL** displays.
14. Press  to save and exit menu to operation mode.
15. Lock the security switch and seal balance if applicable ([Section 5.3 on page 64](#)).

5.2 Perform External Span Calibration Test

A span calibration test is to “check” the current accuracy of the balance. Make sure to use a calibration weight which is equal to the weighing capacity of the balance.





NOTE: Refer to [Table 1-1 on page 7](#) and [Table 1-2 on page 7](#) for model numbers and capacities.

1. Level balance by adjusting the feet, if needed ([Section 2.4 on page 15](#)).
2. Press .  displays.
3. Press  or  to scroll to **6 ADMIN/ADJUST**.
4. Press . **61 SHORT CUT MODE** displays.
5. Press  or  to scroll to **63 MAINTENANCE**.
6. Press . **631 EX CAL** displays.
7. Press  or  to scroll to **632 EX SPAN TEST**.
8. Press . **SPAN TEST STARTED** displays, followed by **ON 0**, then **ON FS** displays.
9. Place full capacity weights on the scale. **ON FS** flashes.
10. Once test is complete, **DATA SAVED** displays, then the instrument error in your selected division size displays.



NOTE: Remove all weight before confirming weight amount, zero calibration is performed immediately following this step.

11. Press . **632 EX SPAN TEST** displays.
12. Press  to save and exit menu to operation mode.
13. Perform a new external span calibration if necessary ([Section 5.1 on page 62](#)).



IMPORTANT: A sealed Legal for Trade balance must be reverified if a new calibration is performed.

5.3 Sealing Instructions

There are two approved methods for sealing a TE Series Balance. Both methods restrict access to the security switch, internal electronics, electrical contacts and Legal for Trade configuration parameters.

5.3.1 Sealing Method 1

Sealing Method 1 involves applying two seal stickers, one over the balance's rear enclosure screw and one over the security switch access hole at the bottom front of the balance.



NOTE: Calibration must be completed before sealing the balance. After locking the security switch, calibration cannot be performed. See [Section 5.0 on page 62](#) to calibration balance.

Method 1 Sealing Instructions

1. Turn off balance and unplug power cable.
2. Remove pre-installed square sticker at the front bottom of the scale.

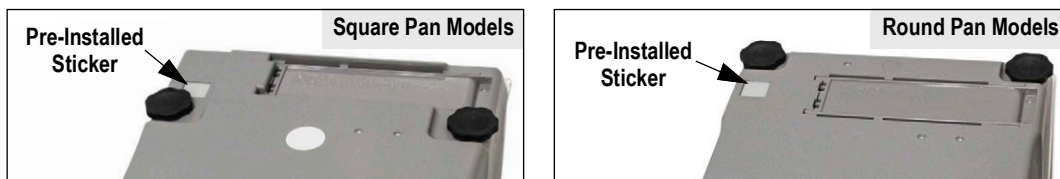


Figure 5-1. Pre-Installed Sticker Location

3. Use a thin non-conductive tool to set the security switch to the locked position (switch initially in unlocked position).



NOTE: Use Sealing Method 2 ([Section 5.3.2 on page 65](#)) if a wire seal is required.

4. Seal balance by placing seal stickers over the rear enclosure screw and the access hole to the security switch.

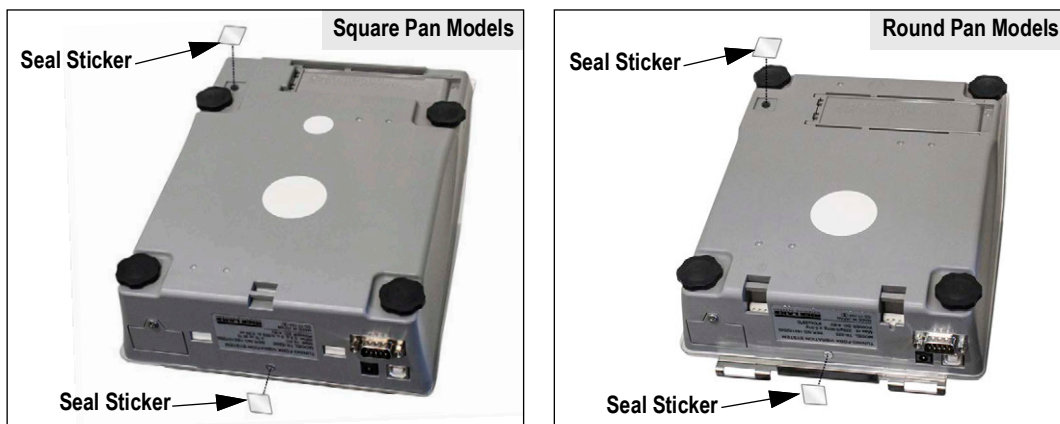


Figure 5-2. Seal Sticker Placement



NOTE: Seal stickers are not included in the box, prepare before inspecting.

5. A weights and measures official or authorized service agent must perform the verification procedure.

5.3.2 Sealing Method 2

Sealing Method 2 involves installing a metal sealing bracket over the balance's rear enclosure screw and applying a seal sticker over the security switch access hole at the bottom front of the balance.



NOTE: Calibration must be completed before sealing the balance. After locking the security switch, calibration cannot be performed. See [Section 5.0 on page 62](#) to calibration balance.

Method 2 Sealing Instructions

1. Turn off balance and unplug power cable.
2. Remove pre-installed square sticker at the front bottom of the scale.

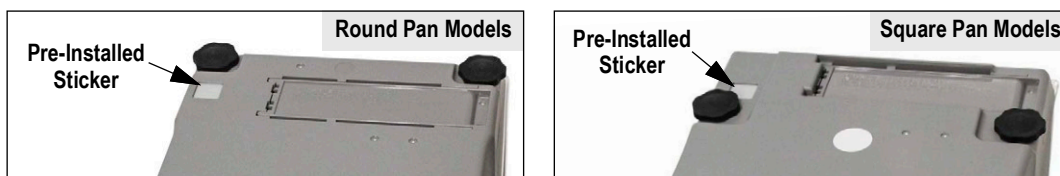


Figure 5-3. Pre-Installed Sticker Location

3. Use a thin non-conductive tool to set the security switch to the locked position (switch initially in unlocked position).
4. Apply a seal sticker over the security switch access hole at the bottom front of the balance.



Figure 5-4. Seal Sticker Placement



NOTE: Seal sticker is not included in the box, prepare before inspecting.

5. For round pan models only, remove screws from top of the base plate with a 5.5 mm hex head screwdriver and remove the base plate.

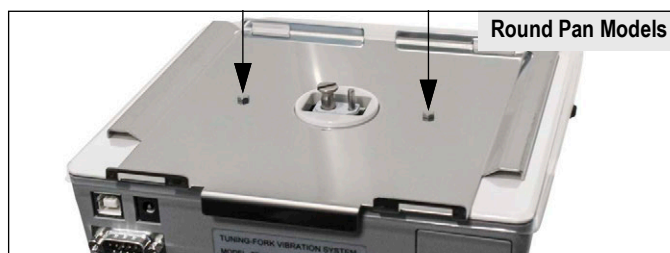


Figure 5-5. Remove Base Plate Screw Locations (Round Pan Models Only)

6. Remove the rear enclosure screw from the back of the balance.



Figure 5-6. Rear Screw Location

7. Secure the metal bracket to the back of the balance with the provided screw included with the bracket.



Figure 5-7. Metal Bracket Placement

8. Position the metal u-bracket inside of the already attached bracket so it covers up the rear screw.
9. Navigate the sealing wire through both sets of bracket holes and seal the wire to secure.



Figure 5-8. Metal U-Bracket and Wire Seal Placement



NOTE: Sealing wire and lead seal are not included in the box, prepare before inspecting.

10. For round pan models only, reposition base plate on top of the balance and secure with a 5.5 mm hex head screwdriver and the previously removed screws.



Figure 5-9. Secure Base Plate Screw Locations (Round Pan Models Only)

11. A weights and measures official or authorized service agent must perform the verification procedure.

6.0 Communications

The balance can be connected to a computer using a compatible third party software program. Connections can be made using RS-232 or USB interfaces. Basic specifications include:

- RS-232C full duplex or USB half duplex
- Asynchronous communication
- RS-232C: EIA-232-D/E or USB: USB 2.0
- Baud Rate: 1200-115200 bps
- Transmission: 1 start bit, non/odd/even number parity, 8 data bit, 1-2 stop bits

6.1 RS-232 Connections

The balance can be equipped with an optional RS-232 feature for communication with printers and computers.

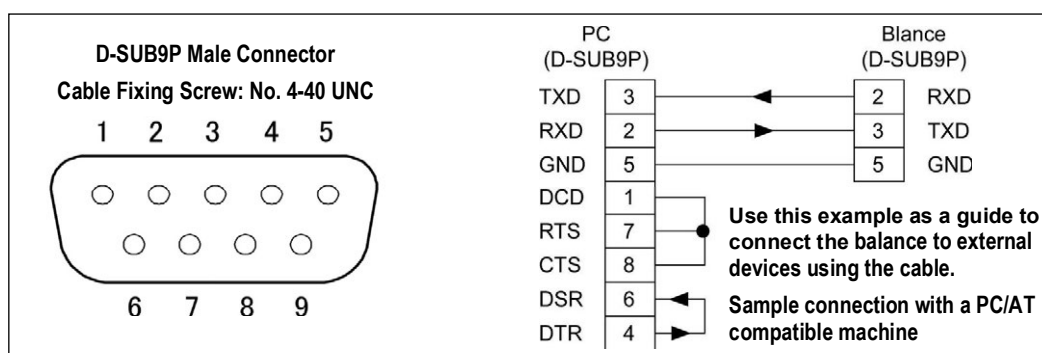


Figure 6-1. RS-232 Connection

Pin Number	Signal	Input/Output	Function
1	-	-	-
2	RXD	Input	Receive Data
3	TXD	Output	Transmit Data
4	DTR	Output	High (when the balance is powered on)
5	GND	-	Signal Grounding
6	-	-	-
7	-	-	-
8	-	-	-
9	External Tare or Zero-Point Adjustment	Input	External Tare Range Setting (NTEP only models); or External Zero-Point Adjustment Setting (NTEP/Measurement Canada models), not available on verified Legal for Trade balances

Table 5-1. Pin Connections



NOTE: The DB-9 connector can set a tare range or adjust the zero point from an external device by connecting a contact or a transistor switch between the pin for externally setting a tare range (pin 9) and the signal ground (pin 5).

Allow at least 400 ms for connection (ON) time (Maximum voltage: 15 V when the balance is turned OFF, sink current: 20 mA when it is turned ON).

6.2 USB Connections

The USB (Type B) connection has the following pin numbering and pin assignments for the balance and connecting cable.

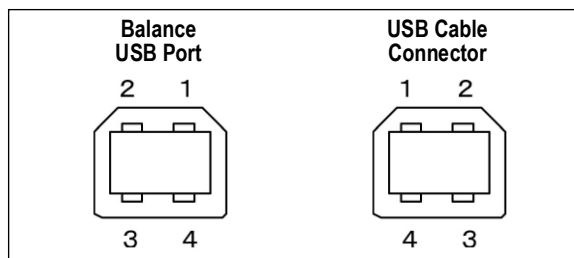


Figure 6-2. USB Connection Pin Numbering

Pin Number	Signal	Function
1	V _{BUS}	Bus power input rating 4.4–5.25 V Cable connector rating 4.75–5.25 V
2	D-	Data signal
3	D+	Data signal
4	GND	Signal grounding

Table 5-2. USB Pin Assignments

6.3 Basic Data Output Format

Date bit: 8 bit, Parity bit/Stop bit: Can be changed.

6-Digit Numeric Format

Consists of 14 characters, including terminators (CR=0xDH/LF=0xAH).

1	2	3	4	5	6	7	8	9	10	11	12	13	14
P1	D1	D2	D3	D4	D5	D6	D7	U1	U2	S1	S2	CR	LF

7-Digit Numeric Format

Consists of 15 characters, including terminators (CR=0xDH/LF=0xAH).

1	2	3	4	5	6	7	8	9	10	11	12	13	14	15
P1	D1	D2	D3	D4	D5	D6	D7	D8	U1	U2	S1	S2	CR	LF

8-Digit Numeric Format

Consists of 16 characters, including terminators (CR=0xDH/LF=0xAH).

1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16
P1	D1	D2	D3	D4	D5	D6	D7	D8	D9	U1	U2	S1	S2	CR	LF

6.3.1 Data Description

Symbol		Code		Description
P1 (one character) indicates the polarity of the data				
+		0x2B		Zero or positive data
-		0x2D		Negative data
D1 to D7/8/D9 (seven or eight or nine characters) stores the numeric data				
0-9		0x30-0x39		0 to 9 (numeric)
.		0x2E		Decimal point (floating)
(SP)		0x20		A space at the top of a numeric value Output to the least significant digit in the absence of a decimal point Unused high-order digit
U1, U2 (two characters) indicates the unit used to show numeric data				
M	G	0x4D	0x47	mg (milligram)
(SP)	G	0x20	0x47	g (gram)
C	T	0x43	0x54	ct (carat)
M	O	0x4D	0x4F	mom (momme)
O	Z	0x4F	0x5A	oz (ounce)
L	B	0x4C	0x42	lb (pound)
O	T	0x4F	0x54	ozt (troy ounce)
D	W	0x44	0x47	dwt (pennyweight)
G	R	0x47	0x52	GN (grain)
T	L	0x54	0x4C	tlH (Hong Kong tael)
T	L	0x54	0x4C	tlS (Singapore, Malaysia tael)
T	L	0x54	0x4C	tlT (Taiwan tael)
t	o	0x74	0x6F	to (tola)
M	S	0x4D	0x53	MSG (mesghal)
B	A	0x42	0x41	BAt (baht)
P	C	0x50	0x43	PCS (parts counting)
(SP)	%	0x20	0x25	% (percentage weighing)
(SP)	#	0x20	0x23	# (multiplied by the coefficient)
(S1) (one character) indicates the judgment result when the limit function is used				
L		0x4C		Shortage (low)
G		0x47		Proper (ok)
H		0x48		Over (high)
(SP)		0x20		No judgment result or data type specified
e		0x65		Net weight
f		0x66		Tare weight; only for NTEP only models
P		0x50		Preset tare weight; only for NTEP only models
T		0x54		Total value (accumulated value)
U		0x55		Unit weight
d		0x64		Gross
S2 (one character) indicates the status				
S		0x53		Data stable
U		0x55		Data unstable
E		0x45		Data error (indicates that data other than S2 is invalid)
(SP)		0x20		No status specified

Table 5-3. Data Description

6.4 CBM Data Output Format

Consists of 26 characters, including terminators (CR=0xDH/LF=0xAH).

(Data bit: 8, Parity: stop, stop bit: can be changed)

1	2	3	4	5	6	7	8	9	10	11	12	13	(SP): space
S1	C1	(SP)	T1	T2	T3	T4	T5	T6	D1	D2	D3	D4	
14	15	16	17	18	19	20	21	22	23	24	25	26	
D5	D6	D7	D8	D9	D10	D11	D12	U1	U2	(SP)	CR	LF	

1	2	3	4	5	6	7	8	9	10	11	12	13	(SP): space
*	*	(SP)	E	R	R	O	R	(SP)	*	*	*	*	
14	15	16	17	18	19	20	21	22	23	24	25	26	
*	*	*	*	*	*	*	*	*	*	(SP)	CR	LF	

6.4.1 Data Description

Symbol	Code	Description
[S1] (1 character) Represents the status.		
(SP)	0x20	Data stable
*	0x2A	Data unstable
[C1] (1 character) Represents the result of comparator function.		
(SP)	0x20	Comparator result: Proper (OK) or No result
H	0x48	Over (HIGH)
L	0x4C	Shortage (LOW)
[T1-T6] (6 characters) Represents the type of the data.		
(SP) (SP) (SP) (SP) (SP) (SP)	0x20 0x20 0x20 0x20 0x20 0x20	Net weight (not tared); only for NTEP only models
N (SP) (SP) (SP) (SP) (SP)	0x4E 0x20 0x20 0x20 0x20 0x20	Net weight (tared); only for NTEP only models
P T (SP) (SP) (SP) (SP)	0x50 0x54 0x20 0x20 0x20 0x20	Preset tare weight; only for NTEP only models
T (SP) (SP) (SP) (SP) (SP)	0x54 0x20 0x20 0x20 0x20 0x20	Tare weight; only for NTEP only models
T O T A L (SP)	0x54 0x4F 0x54 0x41 0x4C 0x20	Total value (Accumulated value); only for NTEP only models
G (SP) (SP) (SP) (SP) (SP)	0x47 0x20 0x20 0x20 0x20 0x20	Gross
U N I T (SP) (SP)	0x55 0x4E 0x49 0x54 0x20 0x20	Unit weight
[D1-D12] (12 characters) Numeric value data is stored.		
+	0x2B	When the data are 0 or positive
-	0x2D	When the data are negative
0-9	0x30 – 0x39	Numeric value (0 – 9)
.	0x2E	Decimal point (floating decimal point)
[0x5B	The number surrounded by '[' and ']' means auxiliary indication
]	0x5D	
(SP)		Spaces fill the top of the data Output to the least significant digit in the absence of a decimal point Unused high-order digit

Table 6-4. Data Description

Symbol		Code		Description	
[U1, U2] (2 characters) Represents the unit of numeric value data.					
m	g	0x6D	0x67	mg	(milligram)
(SP)	g	0x20	0x67	g	(gram)
c	t	0x63	0x74	ct	(carat)
m	o	0x6D	0x6F	mom	(momme)
o	z	0x6F	0x7A	oz	(ounce)
l	b	0x6C	0x62	lb	(pound)
O	T	0x4F	0x54	ozt	(troy ounce)
d	w	0x64	0x77	dwt	(pennyweight)
G	R	0x47	0x52	GN	(grain)
t	l	0x74	0x6C	tlH	(Hong Kong tael)
t	l	0x74	0x6C	tlS	(Singapore, Malaysia tael)
t	l	0x74	0x6C	tlT	(Taiwan tael)
t	o	0x74	0x6F	to	(tola)
M	S	0x4D	0x53	MSG	(mesghal)
B	A	0x42	0x41	BAt	(baht)
P	C	0x50	0x43	PCS	(parts counting)
(SP)	%	0x20	0x25	%	(percentage weighing)
(SP)	#	0x20	0x23	#	(multiplied by coefficient)

Table 6-5. Data Description (continued)

6.5 Input Commands

Input commands can be entered from an external device. [Table 6-6 on page 71](#) displays operation mode input commands.

Commands			
Operation Mode	Zero-point Adjustment Tare Subtraction* Date/time Output	Output Control Comparator Setting Preset Tare Setting* Interval Time Setting	External Contact Input
Weighing	x	x	x
Counting	x	x	x
Percentage	x	x	x
Multiply	x	x	x
Specific gravity	x	-	x
Statistics	x	-	x
Animal	x	-	x
Formulation*	-	-	-

Table 6-6. Transmission Procedure



NOTE: *Formulation Mode, Tare Subtraction, and Preset Tare Setting are only available for NTEP only models.

Select an input command. The balance sends normal completion response or the requested result data.

- The balance transmits an error response if the operation is unsuccessful or if the command is invalid
- In normal display mode, the balance sends a response within one second of receiving the command. A response is sent for tare range (only for NTEP only models), span adjustment or span test commands
- Do not send a command to the balance until a response from the previous command is received from the balance

The balance needs additional response time in some situations:

- The balance waits for stability after receiving a tare (only for NTEP only models) or a zero-point adjustment command if <17 WT STABLE> is <ON>
- If the balance receives a command when setting a function, when it is under span adjustment or if it is busy for other reasons, the command is executed after the current operation is completed

6.5.1 Input Command Composition 1

This is composed of four characters including the terminator (CR=0xDH/LF=0cAH).

C1	C2	CR	LF
----	----	----	----

6.6 Command Formats



IMPORTANT: Do not confuse the alphabetical O for Arabic number 0 (Zero).

C1	C2	Code (C1)	Code (C2)	Description	Response	
					A00/Exx format	ACK/NAK format
T	(SP)	0x54	0x20	Tare Only available for NTEP only models	A00: Normal response	ACK: Normal response
Z	(SP)	0x5a	0x20	Zero-point adjustment	E01: Abnormal response	NAK: Abnormal response
O	0	0x4f	0x30	Stop output		
O	1	0x4f	0x31	Continuous output		
O	2	0x4f	0x32	Continuous output (no output when unstable)		
O	3	0x4f	0x33	Press output key for one-time output		
O	4	0x4f	0x34	Auto output		
O	5	0x4f	0x35	One-time output when stable		
O	6	0x4f	0x36	One-time output when stable		
O	7	0x4f	0x37	Press output key for one-time output when stable		
Commands O0 to O7:						
<ul style="list-style-type: none"> Have the same functions as the output control set by the setting menu Once commands are executed, that state is maintained. The status is reset to the setting menu when the balance is turned on again 						
O	8	0x4f	0x38	One-time output		
O	9	0x4f	0x39	One-time output after stability is obtained		
Command O8 to O9:						
<ul style="list-style-type: none"> Are used to request data from the balance After the command is executed, it returns to O0 						
O	A	0x4f	0x41	Interval function (Output each time the output time has elapsed)		
O	B	0x4f	0x42	Interval function (Output when stable each time the output time has elapsed)		
When OA or OB command is sent, the interval function starts and must be sent again to end the function						

Table 6-7. Zero Point Adjustment/Tare/Output Control Setting Command

C1	C2	Code (C1)	Code (C2)	Description	Response
D	D	0x44	0x44	Date output request	Date data
D	T	0x44	0x54	Time output request	Time data

Table 6-8. Data Output Request and Time Output Request

6.6.1 Input Command Composition 2

This is composed of 15 characters including the terminator (CR=0xDH/LF=0xAH).

1	2	3	4	5	6	7	8	9	10	11	12	13	14	15
C1	C2	,	C3	C3	C3	C3	C3	C3	C3	C3	C3	C3	CR	LF

C3 has a 10 digit maximum (including the polarity +/-, comma and point). Do not include the measuring unit.

Example: Upper limit input 1200.00g: LA,1200.00

Preset tare input 1000.00g: PT,1000.00 (only available for NTEP only models)

Interval time input 12:34:56: IA,12,34,56 (marked off by commas)

The Input command can be entered when Weighing Mode, Percentage Mode, Counting Mode or Multiplied by Coefficient Mode are being used. If input when in an other modes, the outputs an abnormal response.

If the input value is not correct, the balance outputs an abnormal response.

C1	C2	Code (C1)	Code (C2)	Description	C3	Response	
						A00/Exx format	ACK/NAK format
L	A	0x4C	0x41	Lower limit value setting	Numeric value setting	A00: Normal response	ACK: Normal response
L	B	0x4C	0x42	Upper limit value setting	Numeric value setting	E01: Abnormal response	NAK: Abnormal response
L	C	0x4C	0x43	Reference value setting	Numeric value setting	Abnormal response	Abnormal response

Table 6-9. Comparator Setting Command

C1	C2	Code (C1)	Code (C2)	Description	C3	Response	
						A00/Exx format	ACK/NAK format
P	T	0x50	0x54	Preset tare value setting only available for NTEP only models	Numeric value setting	A00: Normal response E01: Abnormal response	ACK: Normal response NAK: Abnormal response

Table 6-10. Preset Tare Value Setting Command



NOTE: When the normal response, the preset tare value is input in 321 PRESET 1 and the balance operates Preset tare. If the input value is 0 at Preset tare setting value command, the preset tare operation is canceled.

C1	C2	Code (C1)	Code (C2)	Description	C3	Response	
						A00/Exx format	A00/Exx format
I	A	0x49	0x41	Interval (output) time setting	Numeric value setting	A00: Normal response E01: Abnormal response	ACK: Normal response NAK: Abnormal response

Table 6-11. Interval (Output) Time Setting Command

6.7 Response

The Response Command Format consists of five characters including terminators.

1	2	3	4	5
A1	A2	A3	CR	LF

Table 6-12. Response Command Format (A00/Exx format)

A1	A2	A3	Code (A1)	Code (A2)	Code (A3)	Description
A	0	0	0x41	0x30	0x30	Normal response
E	0	1	0x45	0x30	0x31	Abnormal response

Table 6-13. Response Command

Response Command Format consists of one character without a terminator.

1
A1

Table 6-14. Response Command Format (ACK/NAK format)

A1	code(A1)	Description
ACK	0x06	Normal response
NAK	0x15	Abnormal response

Table 6-15. Response Command

6.8 External Contact Input

D-sub9 connectors can set a tare range (only for NTEP only models) or adjust the zero-point from an external device by connecting a contact or a transistor switch between the pin for externally setting a tare range (Pin 9, only for NTEP only models) and the signal ground (Pin 5). Allow at least 400 ms for the connection (ON) time.

Maximum voltage: 15 V when the balance is turned off. Sink current: 20 mA when it is turned on.



NOTE: While external contact input is selected, command input is not available.
There is no response command corresponding to external contact input.

6.9 Communication Settings

Set the communication source desired:

- Standard RS232C
- Standard USB Communication
- Extension RS232/Ethernet (Optional)
- Relay Contact Output (Optional)










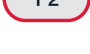


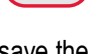




NOTE: *Output condition, 413 CONDITION 1,3,6, cannot be selected.*

433 CONDITION 1,3 and 6 cannot be selected only when Extension RS232C option is connected.










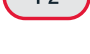
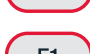


41A STATUS, 42A STATUS and 43A STATUS cannot be selected. The net value status is always appended.



Use the following steps to activate the desired communication source.

1. Press .  displays.
2. Press  or  to select **EXTERNAL I/O**.
3. Press . The current setting displays.
4. Press  or  to select the communication source desired.
5. Press  to enter the menu.
6. Press  or  to display **ACTIVATE**.
7. Press . The current setting begins flashing.
8. Press  or  to select **OFF** or **ON**.
9. Press  to save the setting.
10. Press  to return to operation display.

6.10 Set Communication Parameters

To set the parameters for the communication source:

1. Press .  displays.
2. Press  or  to select **EXTERNAL I/O**.
3. Press . The current setting displays.
4. Press  or  to select the communication source desired.
5. Press  to enter the menu.
6. Press  or  to display the parameter to set.
7. Press . The current parameter setting begins flashing.
8. Press  or  to select desired setting.

9. Press  to save the setting.
10. Repeat process until all parameters have been set.
11. Press  to return to operation display.

See [Section 4.6 on page 54](#) for parameters and settings that are available.

6.10.1 Relay Contact Output (Option)

When using the Relay Contact Output, the **COMPARE** parameter in the **Extension RS232/Ethernet (Optional)** communication source, must be set to:

- 0: AA's per the comparator setting. See [Section 3.15 on page 35](#)
- 1: Output when discrimination result is OK or absent

Set the **EXTERNAL I/O** to **OP RELAY**. See [Section 6.9](#).

6.11 USB Communication and Bus Power Inputs

The balance can communicate using a USB connection or a bus power input.



CAUTION: The internal calibration device cannot be driven with power supplied from the USB.

1. Go to the following website address and click the Downloads tab to download the necessary USB driver.
<https://www.silabs.com/developers/usb-to-uart-bridge-vcp-drivers>



NOTE: If the above website address doesn't work, please visit the Silicon Labs website (<https://www.silabs.com>) and search: "CP210X USB to UART Bridge VCP Drivers".

2. Connect the balance to the PC.
3. Turn on the PC.
4. Configure the communication settings of the PC.
 - Windows® 7 – open the device managers window by going to the Start menu, right clicking the computer, selecting properties and selecting device manager
 - Click the port (COM and LPT) to open the thread and double click the Silicon Labs SP210xUSB to UART Bridge (COM) to open the properties window
 - Go to the Port tab
 - Input the communication setting in accordance with the communication settings of the balance
5. Set the USB power setting of the PC to avoid unexpected shut down of the balance.
 - Windows® 7 – go to the power management tab of the Silicon Labs SP210xUSB to UART Bridge (COM) to open the properties window.
 - Un-check the **Allow the computer to turn off this device to save power** check-box, then press the **OK** key.

7.0 Troubleshooting and Maintenance

This section covers basic troubleshooting and maintenance of the balance.



CAUTION: Observe proper disposal. This balance, including accessories, may not be disposed of in domestic waste in conformance with the specific requirements of the country, county and local jurisdictions. When disposing of the product, contact local authorities and ask for the proper method of disposal.



IMPORTANT: Do not use volatile solvents on the balance.
Unplug the AC adapter from the receptacle when the balance is not going to be used for a long period of time.

7.1 Maintenance Precautions

Dirt or liquids on the weighing pan can cause errors or an unstable weight reading. Clean the balance frequently, ensuring that dust or liquids don't enter into the internal parts of the balance. Use the following instructions to complete general maintenance on the TE balance.

7.1.1 Cleaning—Round Pan Type

1. Remove the draft shield. Refer to [Section 2.5 on page 16](#) to remove the draft shield.

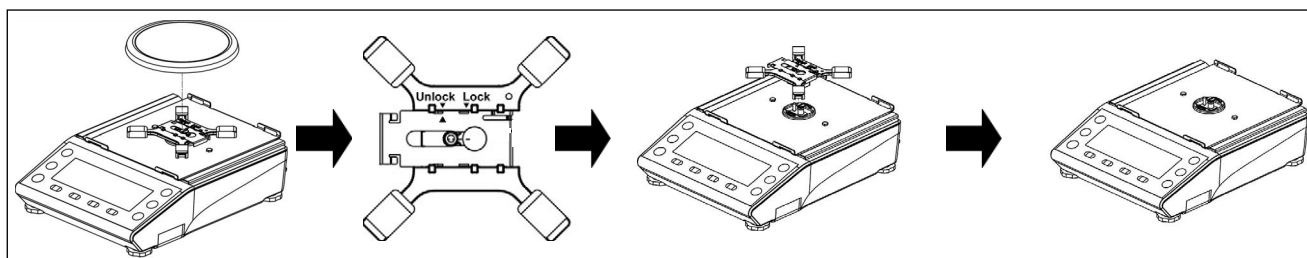


Figure 7-1. Disassemble to Clean

2. Remove the round pan.
3. Move the slider to the unlock side.
4. Remove the pan base.
5. Wipe away dirt with a dry, soft clean cloth. If heavily soiled, remove the weigh pan and the pan base and clean with a clean cloth slightly dampened with a neutral detergent or solvent.

7.1.2 Cleaning—Square Pan Type

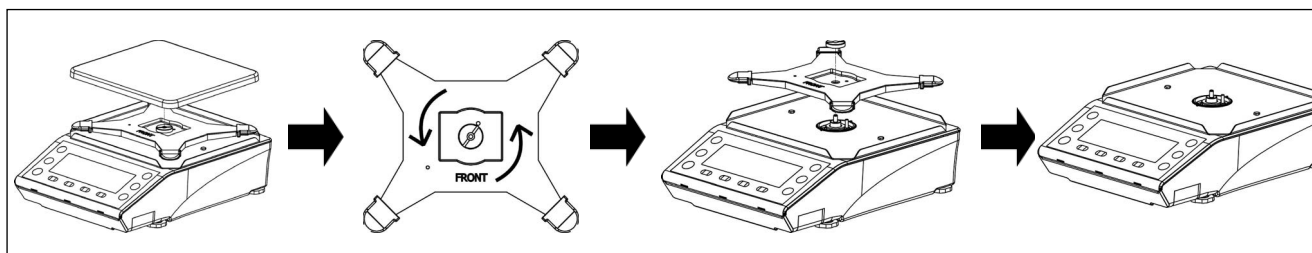


Figure 7-2. Disassemble to Clean

1. Remove the square pan.
2. Remove the pan base.
3. Wipe away dirt with a dry, soft clean cloth. If heavily soiled, remove the weigh pan and the pan base and clean with a clean cloth slightly dampened with a neutral detergent or solvent.

7.2 Error Messages

Error Message/ Error Code	Cause	Solution
OVER ERROR	• Weight exceeds maximum capacity	• Split load into several smaller loads and weigh them • Replace the tare with a lighter one
	• Result exceeds the maximum display digit	• Clear the calculation result, if the calculation results exceed the maximum display digit it will continue to do so
UNDER ERROR	• The load is below the lower limit	• Check the position of the weigh pan and re-seat if necessary • Check for contact with other objects; use the included pan base only
DATA MAX ERROR	• Amount of data exceeds memory	• Clear the data
DISPLAY ERROR	• The result exceeds the maximum display digit	• If calculation results exceed the maximum display digit, it will continue to do so unless something is changed
LOWER ERROR	• The specimen weight/reference weight in Counting/Percentage Mode is below the lower limit	• Ensure the specimen weight/reference weight is higher than the lower limit
ERR001 ERR099	• System error	• Record error code and notify the dealer or Rice Lake Weighing Systems
ERR703	• Operation key pushed when unit is in standby • Hardware issue	• Ensure operation key is not pushed when unit is in standby or starting up from standby • Record error code and notify the dealer or Rice Lake Weighing Systems
ERR705	• Initial zero adjustment error • Initial zero adjustment was not completed during startup due to the unstable load	• Ensure weigh pan is properly seated • Check for contact with other objects • Check for wind or vibration
ERR706	• Load is out of initial zero adjustment range	• Ensure load pan is empty when balance is powered on
ERR709 ERR710 ERR711	• Load is unstable at zero adjustment/tare subtraction • Span adjustment time-out error	• Improper setting of the weighing pan or pan base is suspected • Check for contact with other objects • Check for wind or vibration
ERR717	• Mass of calibration weight differs from designated mass by 1% or more in external span adjustment	• Check the calibration value of the weight and use the proper calibration weight
ERR718	• The mass of the calibration weight is under 50% of the maximum capacity at span adjustment or internal span adjustment weight adjustment by external • calibration weight	• Use a calibration weight which is equal to the maximum capacity
ERR719	• Adjust value by external span adjustment or internal span adjustment is over 1% of maximum capacity	• Execute a 637 REF CAL RESTORE , then execute an internal span adjustment • Check the mass of the weight used for the external span adjustment • Perform a 636 REF CAL
ERR722	• Tare key is pushed during the preset tare operation	• Ensure Tare key is not pushed during preset tare operation
ERR723	• Out of Zero adjustment range (1.5% of maximum capacity)	• Ensure weigh pan is empty while performing a zero adjustment
ERR724	• Out of Tare subtraction range (0 g to maximum capacity)	• Ensure tare weight is within the tare subtraction range
ERR734	• Weight of the sample is out of the importing range at actual value setting method at Percentage Mode (lower limit to maximum capacity)	• Load the sample of which weight is within the importing range
ERR735	• Time-out error of importing the sample weight in the actual value setting method at Percentage Mode	• Improper setting of the weighing pan or pan base is suspected • Check for contact with other object • Check for any wind or vibration
ERR736	• The setting value is out of the setting range at numeric value setting method at Percentage Mode (lower limit to maximum capacity)	• Set the value within the range

Table 7-1. Error Codes

Error Message/ Error Code	Cause	Solution
ERR737	<ul style="list-style-type: none"> Sample weight in the air is out of importing range at Specific Gravity Mode (over 0g to maximum capacity) Sample weight in the water/liquid is out of the importing range at Specific Gravity Mode (0 – maximum capacity to maximum capacity) 	<ul style="list-style-type: none"> Divide the sample so its weight in the air is within the importing range.
ERR738	<ul style="list-style-type: none"> Time-out error in importing the sample weight in the water/liquid at Specific Gravity Mode 	<ul style="list-style-type: none"> Ensure weigh pan is properly seated Check for contact with other object Check for draft or vibration
ERR739	<ul style="list-style-type: none"> Time-out error of importing the sample weight of the value setting method of Preset tare setting 	<ul style="list-style-type: none"> Ensure weigh pan is properly seated Check for contact with other objects Check for any draft or vibration
ERR740	<ul style="list-style-type: none"> The setting value is out of the setting range at numeric value setting method or actual value setting method at Preset tare setting (0g to maximum capacity) 	<ul style="list-style-type: none"> Set the tare within the tare subtraction range
ERR741	<ul style="list-style-type: none"> 631 EX CAL is executed while the external span adjustment function is disabled 	<ul style="list-style-type: none"> Contact the dealer or Rice Lake Weighing Systems
ERR742	<ul style="list-style-type: none"> 633 INT CAL, 634 INT SPAN TEST or 636 REF CAL executed while the balances power is supplied from the USB only 	<ul style="list-style-type: none"> Connect to AC adapter; or insert dry cell batteries and disconnect the USB cable
	<ul style="list-style-type: none"> Internal span adjustment device is not working 	<ul style="list-style-type: none"> Contact the dealer or Rice Lake Weighing Systems
ERR743	<ul style="list-style-type: none"> Battery power is too low to execute 633 INT CAL, 634 INT SPAN TEST or 636 REF CAL 	<ul style="list-style-type: none"> Replace batteries
ERR747	<ul style="list-style-type: none"> When importing a specimen weight in comparator function value setting method, there is a time-out error 	<ul style="list-style-type: none"> Ensure weigh pan is properly seated Check for contact with other objects Check for draft or vibration
ERR748	<ul style="list-style-type: none"> The setting value is out of the setting range at numeric value setting method or actual value setting method at Comparator mode (0 – maximum capacity to maximum capacity) 	<ul style="list-style-type: none"> Set value within range
ERR749	<ul style="list-style-type: none"> When importing a specimen weight in adding function value setting method, there is a time-out error 	<ul style="list-style-type: none"> Ensure weigh pan is properly seated Check for contact with other objects Check for draft or vibration
ERR750	<ul style="list-style-type: none"> Weight of added sample is out of range (0 – maximum capacity) 	<ul style="list-style-type: none"> Choose sample within importing range
	<ul style="list-style-type: none"> The total value exceeds maximum display digit 	<ul style="list-style-type: none"> Clear total value
ERR751	<ul style="list-style-type: none"> Weight of sample is lighter than the minimum interval of the balance in Counting Mode 	<ul style="list-style-type: none"> Choose a sample which unit weight is more than minimum interval of the balance
ERR752	<ul style="list-style-type: none"> Weight of sample is 0g or under in Counting Mode 	<ul style="list-style-type: none"> Choose sample which unit weight is more than the minimum interval of the balance Counting Mode cannot operate subtractive counting
ERR753	<ul style="list-style-type: none"> Time-out error of importing the unit weight at Counting Mode 	<ul style="list-style-type: none"> Ensure weigh pan is properly seated Check for contact with other objects Check for draft or vibration
ERR754	<ul style="list-style-type: none"> Deleted the latest data then executed deleting operation of the second latest data at Statistics Mode 	<ul style="list-style-type: none"> Only the latest data can be deleted Select <ALL> to delete all the other data
ERR755	<ul style="list-style-type: none"> Time-out error of importing the sample weight at Statistics/Formulation Mode 	<ul style="list-style-type: none"> Ensure weigh pan is properly seated Check for contact with other objects Check for draft or vibration
ERR756	<ul style="list-style-type: none"> Weight of the sample is out of the importing range at Statistics/Formulation Mode (0g to maximum capacity) 	<ul style="list-style-type: none"> Choose sample which weighs within the importing range
ERR757	<ul style="list-style-type: none"> Bluetooth connection error 	<ul style="list-style-type: none"> Disconnect and reconnect the Bluetooth communication

Table 7-1. Error Codes (Continued)

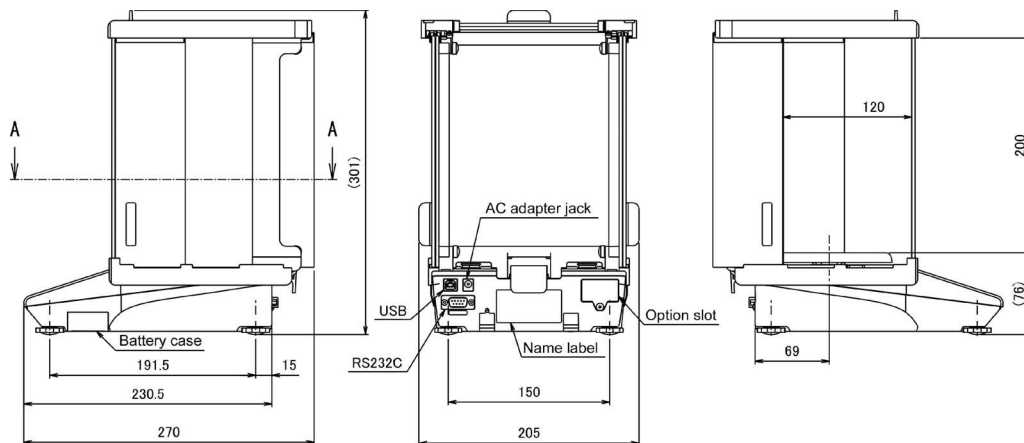
Error Message/ Error Code	Cause	Solution
ERR758	<ul style="list-style-type: none"> Bluetooth hardware error 	<ul style="list-style-type: none"> Contact dealer or Rice Lake Weighing Systems
ERR760	<ul style="list-style-type: none"> Adding operation is executed while the Adding function is disabled 	<ul style="list-style-type: none"> Set 141 ACTIVATE ON then execute the adding operation
ERR761	<ul style="list-style-type: none"> An error occurred at 636 REF CAL 	<ul style="list-style-type: none"> Re-execute 636 REF CAL
ERR763	<ul style="list-style-type: none"> Calculation error in specific gravity of the sample in Specific Gravity Mode 	<ul style="list-style-type: none"> Re-execute the specific gravity function
ERR764	<ul style="list-style-type: none"> External weight used for 631 EX CAL is different from selected weight range at SELECT WEIGHT 	<ul style="list-style-type: none"> Use a weight that is within the selected range

Table 7-1. Error Codes (Continued)

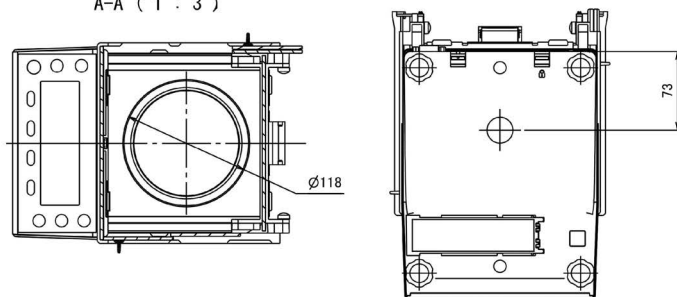
8.0 Specifications

8.1 Dimensions

TE-223
TE-623
TE-322NC



A-A (1 : 3)



TE-3202
TE-6202
TE-15001
TE-1501NC
TE-8200NC

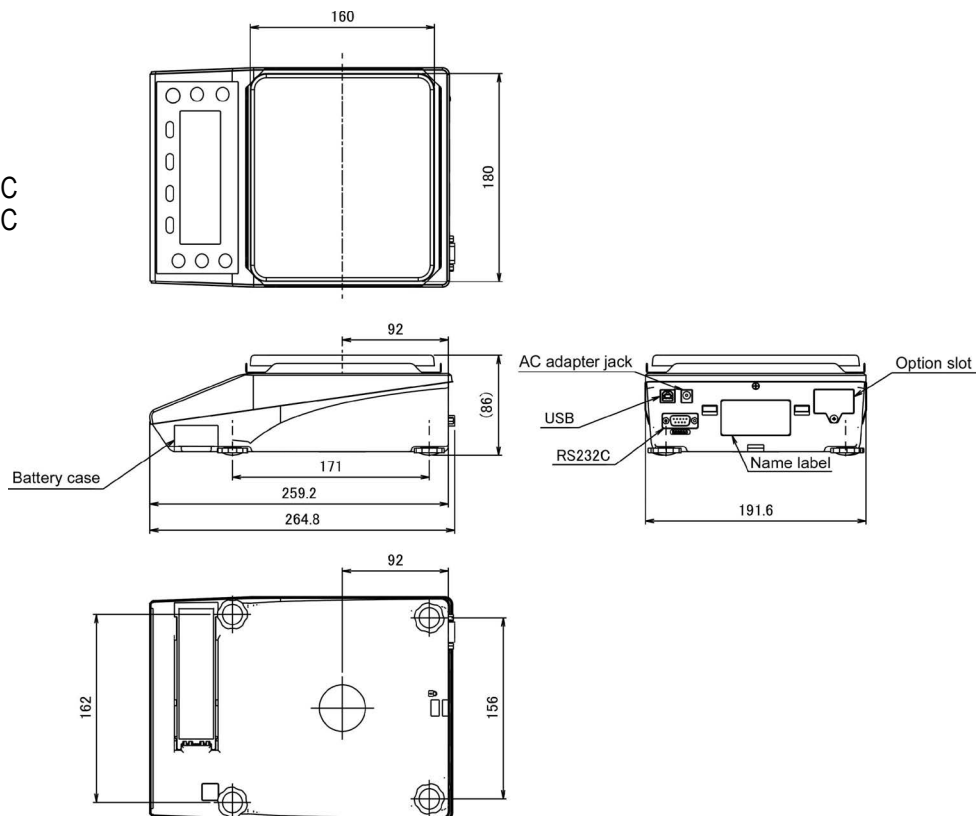


Figure 8-1. TE Model Dimensions

8.2 Model Specifications

Model	Max (g)	e (g)	d (g)	Weighing Range (g)	Accuracy Class	Draft Shield	Span Adjustment
TE-223	220	0.01	0.001	0 - 220.09	II	Yes	External
TE-623	620	0.01	0.001	0 - 620.09			
TE-3202	3200	0.1	0.01	0 - 3200.9			
TE-6202	6200	0.1	0.01	0 - 6200.9		No	
TE-15001	15000	1	0.1	0 - 15009			
TE-322NC	320	0.01	0.01	0 - 320.09		Yes	
TE-1501NC	1500	0.1	0.1	0 - 1500.9			
TE-8200NC	8200	1	1	0 - 8209		No	

Table 8-1. Model Specifications



NOTE: For NTEP only models, $e \neq d$. The balance's readability (d) is 10x the NTEP verified readability (e).
For NTEP/Measurement Canada models, $e = d$. The balance's readability (d) is equal to the verified readability (e).

8.3 Basic Communication Specifications

Communication Method

RS-232C Full-duplex communication method

USB Half-duplex communication method

Synchronization Method

Asynchronous communication method

Electrical Specification

RS-232C EIA-232-D/E

USB USB2.0

Baud rate

1200/2400/4800/9600/19200/38400/57600/115200 bps

Transmission Code Composition

Start bit. 1 bit

Parity bit. None/Odd number/Even number

Data bit. 8 bit

Stop bit. 1 bit/2 bit

Radio Compliance

FCC ID. PVH0946

IC 5325A-0946

8.4 Certifications and Approvals



NTEP

CoC Number: 17-104A1

All TE models are NTEP Certified



Measurement Canada

Approval Number: AM-6169

Only TE-322, TE-1501NC, TE-8200NC
are Measurement Canada Approved

8.5 Functional Specifications

Item	Description
Weighing System	Tuning fork vibration method
Weighing Modes	Weighing / Counting / Percentage / Multiplied by Coefficient** / Animal** / Specific Gravity / Statistics** / Formulation* modes
Functions	Functions related to the operation — Unit Setting / Comparator / Adding** / Tare-Subtraction Reminder* / Zero-Point Adjustment Reminder* / Stability Waiting** / Bar Graph / Back Light / Auto Power Off / Simple SCS
	Function related to the performance** — Stability Discrimination Width** / Response Speed** / Automatic Zero Tracking**
	User information settings — Preset Tare* / Preset Tare Weight* / Compare Weight / Compare Percentage / Compare Counting / Multiplied by Coefficient Comparator**
	Functions related to the lock — Total Lock Release / Key Lock / Menu Lock
	Controlling and adjustment functions — Shortcut / Free Key / Span Adjustment with External Weight** / Span Test with External Weight / Balance ID / Password / Output Language (English, German, Spanish, French, Japanese) / Date Setting / Time Setting / Designation of Minimum Indication / Readability Setting* / Span Adjustment at Power On* / Direct Start
Display	LCD with back light 7-segment: maximum 8-digit / segment height up to 16.5 mm 16-segment: maximum 20-digit / segment height up to 8.5 mm Bar graph: 40-step (NTEP only models) / 30-step (NTEP / Measurement Canada models)
Tare Range Setting	Weight subtraction with the tare key*
Automatic Zero Tracking	Enabled (can be disabled in settings)**
Display When Overloaded	When the indication limit is exceeded, <OVER ERROR> is indicated
Output	RS-232C compliant output is equipped as standard (D-sub9P male connector) USB (type B connector)
Span Adjustment	External span adjustment and calibration
Counting Mode Minimum Unit Weight	TE-223, TE-623: 0.001 g (NTEP 0.01 g) TE-322NC 0.01 g TE-3202, TE-6202: 0.01 g (NTEP 0.1 g) TE-1501NC 0.1 g TE-15001: 0.1 g (NTEP 1 g) TE-8200NC 1 g
Percentage Mode Weight Limit	TE-223 – TE-623: 0.1 g TE-322NC 1 g TE-3202 – TE-6202: 1 g TE-1501NC 10 g TE-15001: 10 g TE-8200NC 100 g
Power	Dedicated AC adapter (100-240 VAC / 50-60 Hz) Dry cell batteries USB bus power: connected with PC in which the driver is installed
Ratings	AC adapter: 4-6 VDC (0.3 A) Battery box (4 AA): 4-6 VDC (0.3 A) USB: 5 VDC (0.3 A)
Dimensions of Weighing Pan	TE-223 – TE-623: 118 mm diameter TE-3202 – TE-15001: 160 x 180 mm
Balance Weight (net)	TE-223 – TE-623: 2.6 kg TE-1502 – TE-15001: 2.7 kg
Operating Conditions	Temperature: 5-35 °C Humidity: 85% or lower (no condensation) Pollution degree: 2 Altitude: 2000 m or less above sea level Location of use: Indoor
Option	Extension RS-232C, Relay Contact, Ethernet
*Only available in NTEP only models; **Not available in verified Legal for Trade NTEP / Measurement Canada models	

Table 8-2. Functional Specifications



© Rice Lake Weighing Systems Content subject to change without notice.

230 W. Coleman St. • Rice Lake, WI 54868 • USA USA: 800-472-6703 • International: +1-715-234-9171